



TOTAL
Solution

INHEALTH
TECHNOLOGIES®

Blom-Singer®

2025 PRODUCT CATALOG

CONTENTS

| | |
|--|-----------|
| Pre-Operative Education | 2 |
| "What to Expect" Brochure | 2 |
| Educational Videos | 2 |
| Surgical Solutions | 3 |
| Tracheoesophageal Puncture Set (Sterile) | 3 |
| Indwelling Voice Prosthesis (Sterile) | 3 |
| Adjustable Bi-Flanged Fistula Prosthesis (Sterile) | 3 |
| Post-Operative Care | 4 |
| Laryngectomy Tubes | 4 |
| Marpac Tracheostomy Tube Holder/Collar | 4 |
| StomaSoft Cleaning Swabs | 4 |
| Electrolarynx | 5 |
| Post-Operative Care Kits (For Total Laryngectomy) | 6 |
| Next Steps | 8 |
| Voice Prosthesis (Clinician Placed) | 8 |
| Special Order Voice Prosthesis (Clinician Placed) | 10 |
| Voice Prosthesis (Patient Changeable) | 12 |
| Voice Prosthesis Accessories | 13 |
| Heat & Moisture Exchange (HME) | 14 |
| Adhesive Housings | 16 |
| Tracheostoma Attachments | 16 |
| And Beyond | 18 |
| Hands-Free Speech Valves | 18 |
| Skin Care | 18 |
| Stoma Protection and Accessories | 19 |
| Professional Products | 21 |
| Voice Prosthesis Management | 22 |
| Oto-Rhino Foreign Body Remover | 24 |
| Nasal Septal Perforation Prosthesis | 26 |
| Ordering Information | 27 |
| New Patient Enrollment Instructions | 27 |
| Free Educational Resources | 27 |

InHealth Technologies® offers a TOTAL Solution for your total laryngectomy from pre- to post-op and beyond

TOTAL Solution serves patients at all phases of voice restoration, from products that deliver solutions for pre- and post-operative care and beyond. This continuum of care is enhanced by InHealth Technologies' customer-focused team, including leading clinical specialists, experienced patient consultants and dedicated customer service members.

Cover Image:

Left: Janice
Year of TL: 2014
Wearing Accufit housing with
SpeakFree HME Hands-Free Valve

Right: Mackayla
Year of TL: 2021
Wearing Hydrofit housing
and Day&Night HME

TOTAL Solution

Pre-Operative Education

- "What to Expect" Brochure
- Educational Videos

PHASE

1

PHASE

2

Surgical Solutions

- Tracheoesophageal Puncture Set
- Sterile Option Voice Prosthesis
- Fistula Prosthesis

PHASE

3

Post-Operative Care

- Post-Operative Care Kits
- Laryngectomy Tubes
- Electrolarynx

PHASE

4

Next Steps

- Next Steps Kit
- Voice Prosthesis
- HMEs
- Adhesive Housings

PHASE

5

And Beyond



- Hands-Free Speech Valve
- Skin Care

● PRE-OPERATIVE EDUCATION

“What to Expect” Brochure

Available online and free to download at:

<https://inhealth.com/educational-materials>

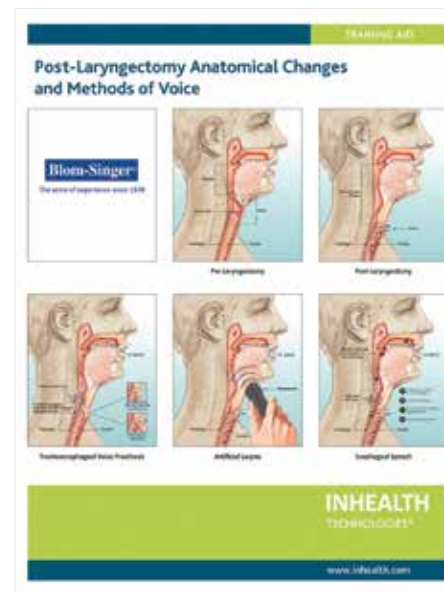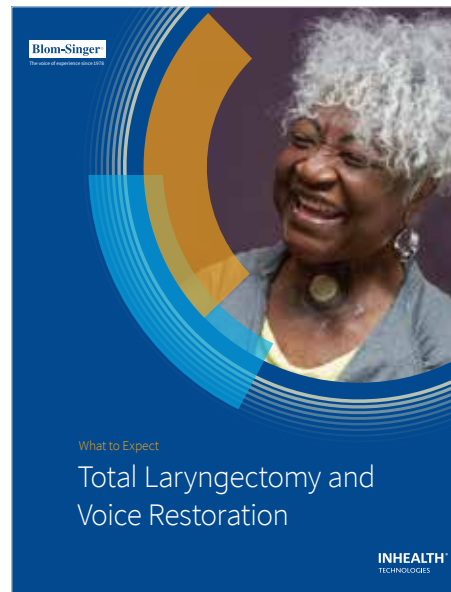

Educational Videos

Available at InHealth Technologies®

<https://inhealth.com/videos>



● SURGICAL SOLUTIONS

Blom-Singer® Tracheoesophageal Puncture Set (Sterile)

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

Designed for secondary placement of a tracheoesophageal prosthesis. Includes a disposable 16 gauge curved introducer needle with a 14 gauge Teflon® sheath, a 26" (66cm) silicone catheter with a tapered reinforced tip attached to a 24" (61cm) Teflon-coated stainless steel leader cable, and a fenestrated surgical drape with adhesive around the aperture.

| Ref | Description |
|---------|--|
| BE 6055 | Tracheoesophageal Puncture Set (Sterile) |

Blom-Singer® Classic™ Indwelling Voice Prosthesis (Sterile)

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

Designed for primary TE puncture and placement procedures. Includes insertion and cleaning accessories. Clinician placed.


| Length | 16 Fr | 20 Fr | Length | 16 Fr | 20 Fr |
|--------|------------|------------|--------|------------|------------|
| 8mm | IN 1608-SO | IN 2008-SO | 12mm | IN 1612-SO | IN 2012-SO |
| 10mm | IN 1610-SO | IN 2010-SO | 14mm | IN 1614-SO | IN 2014-SO |

Blom-Singer® Adjustable Bi-Flanged Fistula Prosthesis (Sterile)

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

Designed by Eric D. Blom, Ph.D to control leakage of saliva, food, and drink around a hypopharyngeal fistula. The prosthesis is adjustable in length and comes in two flange sizes, 25mm and 38mm, ensuring a proper fit.

| Ref | Description |
|---------|--|
| FM 1003 | 25mm (Small) Adjustable Bi-Flanged Fistula Prosthesis, 1 ea |
| FM 1004 | 38mm (Medium) Adjustable Bi-Flanged Fistula Prosthesis, 1 ea |

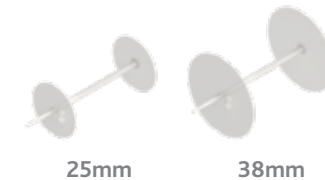

● POST-OPERATIVE CARE

Blom-Singer® StomaSoft® Laryngectomy Tubes

Blom-Singer StomaSoft Laryngectomy Tubes are flexible, soft, and compatible with HME and other accessories from InHealth Technologies®. They feature reinforced strap holders and are available with or without fenestrations.

| Size/Length | I.D./O.D. | Non-Fenestrated | Fenestrated | Rx |
|-------------|---------------|-----------------|-------------|----|
| 8/36 | 9.5mm/12.0mm | BE 6698 | BE 6698F | ■ |
| 8/55 | 9.5mm/12.0mm | BE 6699 | BE 6699F | ■ |
| 9/36 | 10.5mm/13.5mm | BE 6700 | BE 6700F | ■ |
| 9/55 | 10.5mm/13.5mm | BE 6701 | BE 6701F | ■ |
| 10/36 | 12.0mm/15.0mm | BE 6702 | BE 6702F | ■ |
| 10/55 | 12.0mm/15.0mm | BE 6703 | BE 6703F | ■ |
| 12/36 | 13.5mm/17.0mm | BE 6704 | BE 6704F | ■ |
| 12/55 | 13.5mm/17.0mm | BE 6705 | BE 6705F | ■ |



Blom-Singer® Laryngectomy Tubes (Sterile)

Blom-Singer Laryngectomy Tubes (Sterile) are durable, lightweight silicone tubes that are used to keep the tracheostoma open after laryngectomy. The tubes do not hold an HME.

| Size/Length | I.D./O.D. | Non-Fenestrated (Sterile) | Rx |
|-------------|---------------|---------------------------|----|
| 8/36 | 9.5mm/12.0mm | BE 6398 | ■ |
| 8/55 | 9.5mm/12.0mm | BE 6399 | ■ |
| 9/36 | 10.5mm/13.5mm | BE 6400 | ■ |
| 9/55 | 10.5mm/13.5mm | BE 6401 | ■ |
| 10/36 | 12.0mm/15.0mm | BE 6402 | ■ |
| 10/55 | 12.0mm/15.0mm | BE 6403 | ■ |
| 12/36 | 13.5mm/17.0mm | BE 6404 | ■ |
| 12/55 | 13.5mm/17.0mm | BE 6405 | ■ |



Marpac Tracheostomy Tube Holder/Collar™

Marpac Tracheostomy Tube Holder/Collar (for laryngectomy patients) is a single use, disposable Comfortable Collar™ used to secure the tube. Latex free. Adult—one size fits all.

| Ref | Description |
|-----|---------------------------------|
| 204 | Tracheostomy Tube Holder/Collar |



StomaSoft Cleaning Swabs

The StomaSoft Cleaning Swabs are intended for cleaning tubes after removal from the patient. Available in packs of 10.

| Ref | Description |
|---------|-----------------------|
| TSWB-10 | Foam Swab, pack of 10 |



Blom-Singer® Electrolarynx EL 2000

An electrolarynx is a speaking aid that produces sound vibrations that can be conducted through the throat and mouth to speak after a laryngectomy.

InHealth Technologies® and Eric D. Blom, Ph.D are pleased to offer a high quality, cost effective electrolarynx for laryngectomees.

| Ref | Description |
|---------|---|
| EL 2000 | Blom-Singer® Electrolarynx Digital Speech Aid |
| EL 2002 | Oral Adapter |

Oral Adaptor EL 2002



Electrolarynx EL 2000



Blom-Singer® Post-Operative Care Kits (For Total Laryngectomy)

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

The Blom-Singer® Post-Operative Care Kits were designed with valued input from clinicians and patients. These kits contain supplies to assist patients in their transition from the hospital to home after a total laryngectomy. They are available in several configurations to meet select needs. All configurations include your choice of StomaSoft Laryngectomy Tube and HMEs as listed below.

The Post-Operative Care Kit is housed in a hard, clear storage case, which also includes information on how to obtain our Next Steps kit. The Next Steps kit includes additional sample product and educational items.

All kit configurations contain the following patient care supplies:

- 7 ea. 204 Marpac Tracheostomy Tube Holder/Collar™
- BE 6048 Blom-Singer® Shower Guard
- Drawstring Bag
- HydroFit Adhesive Housings, 6 pack
- Mirror with Stand
- Halcyon® LED Flashlight
- Next Steps Kit - Request Flyer



| REF | Description |
|----------------|---|
| PK 2100 | |
| BE 6701 | Blom-Singer® StomaSoft® Laryngectomy Tube 9/55 Non-fenestrated, 2 ea |
| BE 1091 | Day&Night HME® Cartridge ClassicFlow Heat and Moisture Exchange, 1 ea |
| PK 2101 | |
| BE 6703 | Blom-Singer® StomaSoft® Laryngectomy Tube 10/55 Non-fenestrated, 2 ea |
| BE 1091 | Day&Night HME® Cartridge ClassicFlow Heat and Moisture Exchange, 1 ea |
| PK 2102 | |
| BE 6701 | StomaSoft Laryngectomy Tube 9/55 Non-fenestrated, 2 ea |
| BE 1092 | Day&Night HME® Cartridge EasyFlow Heat and Moisture Exchange, 1 ea |
| PK 2103 | |
| BE 6703 | Blom-Singer® StomaSoft® Laryngectomy Tube 10/55 Non-fenestrated, 2 ea |
| BE 1092 | Day&Night HME® Cartridge EasyFlow Heat and Moisture Exchange, 1 ea |
| PK 2104 | |
| EL 2000 | Blom-Singer® Electrolarynx, 1 ea |
| BE 6703 | Blom-Singer® StomaSoft® Laryngectomy Tube 10/55 Non-fenestrated, 2 ea |
| BE 1091 | Day&Night HME® Cartridge ClassicFlow Heat and Moisture Exchange, 1 ea |
| PK 2106 | |
| BE 6705 | Blom-Singer® StomaSoft® Laryngectomy Tube 12/55 Non-fenestrated, 2 ea |
| BE 1091 | Day&Night HME® Cartridge ClassicFlow Heat and Moisture Exchange, 1 ea |
| PK 2107 | |
| BE 6705 | Blom-Singer® StomaSoft® Laryngectomy Tube 12/55 Non-fenestrated, 2 ea |
| BE 1092 | Day&Night HME® Cartridge EasyFlow Heat and Moisture Exchange, 1 ea |

Mike
Year of TL: 2015
Wearing StomaSoft
laryngectomy tube and
Day&Night HME



● NEXT STEPS

VOICE PROSTHESES (CLINICIAN PLACED)

Blom-Singer® Classic™ Indwelling Voice Prostheses (Sterile)

Includes sterile prostheses for primary TE puncture and placement procedures. Includes insertion and cleaning accessories. Clinician placed.

| Length | 16 Fr | 20 Fr | Rx | Length | 16 Fr | 20 Fr | Rx |
|--------|------------|------------|----|--------|------------|------------|----|
| 8mm | IN 1608-SO | IN 2008-SO | ■ | 12mm | IN 1612-SO | IN 2012-SO | ■ |
| 10mm | IN 1610-SO | IN 2010-SO | ■ | 14mm | IN 1614-SO | IN 2014-SO | ■ |



Blom-Singer® Classic™ Indwelling Voice Prostheses

Includes nonsterile voice prosthesis for secondary or replacement TEP procedures, with insertion/cleaning accessories. Clinician placed.

| Length | 16 Fr | 20 Fr | Rx | Length | 16 Fr | 20 Fr | Rx |
|--------|------------|------------|----|--------|------------|------------|----|
| 4mm | IN 1604-NS | IN 2004-NS | ■ | 14mm | IN 1614-NS | IN 2014-NS | ■ |
| 6mm | IN 1606-NS | IN 2006-NS | ■ | 16mm | IN 1616-NS | IN 2016-NS | ■ |
| 8mm | IN 1608-NS | IN 2008-NS | ■ | 18mm | IN 1618-NS | IN 2018-NS | ■ |
| 10mm | IN 1610-NS | IN 2010-NS | ■ | 20mm | IN 1620-NS | IN 2020-NS | ■ |
| 12mm | IN 1612-NS | IN 2012-NS | ■ | | | | |



Blom-Singer® Dual Valve™ Indwelling Voice Prostheses

The Blom-Singer Dual Valve Indwelling Voice Prosthesis has two valves embedded with silver oxide as a material preservative. When the primary valve fails, the secondary speech valve is intended to remain viable to help prevent leakage of fluids. The Blom-Singer® Dual Valve Voice Prosthesis is available with standard or large tracheal/esophageal flanges. Includes nonsterile voice prosthesis for secondary or replacement TEP procedures, with insertion/cleaning accessories. Clinician placed.

| Length | 20 Fr | 20 Fr Large Flange | Rx | Length | 20 Fr | 20 Fr Large Flange | Rx |
|--------|------------|--------------------|----|--------|------------|--------------------|----|
| 6mm | DV 2006-NS | DV 2006-LF | ■ | 12mm | DV 2012-NS | DV 2012-LF | ■ |
| 8mm | DV 2008-NS | DV 2008-LF | ■ | 14mm | DV 2014-NS | DV 2014-LF | ■ |
| 10mm | DV 2010-NS | DV 2010-LF | ■ | | | | |



Blom-Singer® Advantage® Indwelling Voice Prostheses

The Blom-Singer Advantage Indwelling Voice Prosthesis is embedded with silver oxide as a material preservative.* Includes nonsterile voice prosthesis for secondary or replacement TEP procedures, with insertion/cleaning accessories.

HARD VALVE ASSEMBLY

| Length | 20 Fr | Rx |
|--------|------------|----|
| 4mm | AD 2004-NS | ■ |
| 6mm | AD 2006-NS | ■ |
| 8mm | AD 2008-NS | ■ |
| 10mm | AD 2010-NS | ■ |
| 12mm | AD 2012-NS | ■ |
| 14mm | AD 2014-NS | ■ |

SOFT VALVE ASSEMBLY

| Length | 16 Fr | 20 Fr | Rx |
|--------|------------|------------|----|
| 6mm | AS 1606-NS | AS 2006-NS | ■ |
| 8mm | AS 1608-NS | AS 2008-NS | ■ |
| 10mm | AS 1610-NS | AS 2010-NS | ■ |
| 12mm | AS 1612-NS | AS 2012-NS | ■ |
| 14mm | AS 1614-NS | AS 2014-NS | ■ |



*Silver oxide, a known anti-microbial substance, has been incorporated into the silicone flap valve of the Blom-Singer Advantage and Dual Valve Indwelling Voice Prostheses as a material preservative.

Ksenia
Year of TL: 2019
Wearing Accufit
housing and
Day&Night HME



SPECIAL ORDER VOICE PROSTHESES (CLINICIAN PLACED)


—DEVELOPED BY ERIC D. BLOM, PH.D

Unique problems call for unique solutions. InHealth Technologies® offers several voice prostheses to address those needs.

Special Length

Available in one millimeter incremental lengths for a more precise TE puncture fit.


| | 5mm | 7mm | 9mm | Rx |
|-------|------------|------------|------------|----|
| 16 Fr | IN 1605-SL | IN 1607-SL | IN 1609-SL | ■ |
| 20 Fr | IN 2005-SL | IN 2007-SL | IN 2009-SL | ■ |

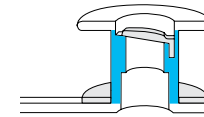

Increased Resistance

Designed for patients who experience intermittent valve opening with fluctuations in esophageal pressures, resulting in air ingestion and/or leakage.


| | 6mm | 8mm | 10mm | 12mm | Rx |
|-------|------------|------------|------------|------------|----|
| 16 Fr | IN 1606-IR | IN 1608-IR | IN 1610-IR | IN 1612-IR | ■ |
| 20 Fr | IN 2006-IR | IN 2008-IR | IN 2010-IR | IN 2012-IR | ■ |

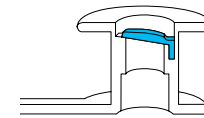

Large Esophageal Flange


Provides a seal to help manage leakage around a voice prosthesis. Also, it enables downsizing from a wider diameter prosthesis to a smaller diameter. Effectiveness depends on an optimal prosthesis length fit that achieves a snug flange-to-flange seal against the wall surface(s).

| | 4mm | 6mm | 8mm | 10mm | 12mm | 14mm | Rx |
|-------|-------------|-------------|-------------|-------------|-------------|-------------|----|
| 16 Fr | IN 1604-LEF | IN 1606-LEF | IN 1608-LEF | IN 1610-LEF | IN 1612-LEF | IN 1614-LEF | ■ |
| 20 Fr | IN 2004-LEF | IN 2006-LEF | IN 2008-LEF | IN 2010-LEF | IN 2012-LEF | IN 2014-LEF | ■ |

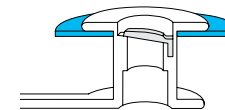

LARGE ESOPHAGEAL FLANGE (SPECIAL LENGTH)

| | 5mm | 7mm | Rx |
|-------|--------------|--------------|----|
| 16 Fr | IN 1605-LESL | IN 1607-LESL | ■ |
| 20 Fr | IN 2005-LESL | IN 2007-LESL | ■ |

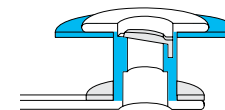

Large Esophageal and Tracheal Flanges

Provides a double seal in patients with enlarged puncture tracts and tear-drop or irregularly shaped punctures. Effectiveness depends on an optimal prosthesis length fit that achieves a snug flange-to-flange seal against the wall surface(s).

| | 4mm | 6mm | 8mm | 10mm | 12mm | 14mm | Rx |
|-------|------------|------------|------------|------------|------------|------------|----|
| 16 Fr | IN 1604-LF | IN 1606-LF | IN 1608-LF | IN 1610-LF | IN 1612-LF | IN 1614-LF | ■ |
| 20 Fr | IN 2004-LF | IN 2006-LF | IN 2008-LF | IN 2010-LF | IN 2012-LF | IN 2014-LF | ■ |




LARGE ESOPHAGEAL FLANGE (INCREASED RESISTANCE)

| | 4mm | 6mm | 8mm | Rx |
|-------|--------------|--------------|--------------|----|
| 16 Fr | IN 1604-LEIR | IN 1606-LEIR | IN 1608-LEIR | ■ |
| 20 Fr | IN 2004-LEIR | IN 2006-LEIR | IN 2008-LEIR | ■ |

TEP Occluders

The Blom-Singer® TEP (Tracheoesophageal Puncture) Occluder is a non-valved silicone device designed for placement in, and maintenance of, the TEP following total laryngectomy. Not a voice prosthesis.

| | 4mm | 6mm | 8mm | 10mm | 12mm | 14mm | Rx |
|-------|--------------|--------------|--------------|--------------|--------------|--------------|----|
| 16 Fr | IN 1604-LETO | IN 1606-LETO | IN 1608-LETO | IN 1610-LETO | IN 1612-LETO | IN 1614-LETO | ■ |
| 20 Fr | IN 2004-LETO | IN 2006-LETO | IN 2008-LETO | IN 2010-LETO | IN 2012-LETO | IN 2014-LETO | ■ |

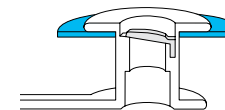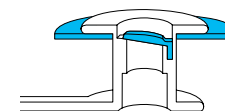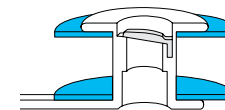

VOICE PROSTHESES (PATIENT CHANGEABLE)

Blom-Singer® Duckbill Voice Prostheses

The Blom-Singer® Duckbill is a patient changeable, slit valve prosthesis.

| Length | 16 Fr | Rx | Length | 16 Fr | Rx |
|--------|-----------|----|--------|-----------|----|
| 6mm | DB 16-006 | ■ | 12mm | DB 16-012 | ■ |
| 8mm | DB 16-008 | ■ | 14mm | BE 6001 | ■ |
| 10mm | BE 6000 | ■ | 18mm | BE 6002 | ■ |

Blom-Singer® Low Pressure Voice Prostheses

The Blom-Singer® Low Pressure is a patient changeable, one-way valve prosthesis.

| Length | 16 Fr | 20 Fr | Rx | Length | 16 Fr | 20 Fr | Rx |
|--------|-----------|-----------|----|--------|---------|---------|----|
| 6mm | BE 6009 | BE 2014 | ■ | 18mm | BE 6012 | BE 2026 | ■ |
| 8mm | LP 16-008 | LP 20-008 | ■ | 22mm | BE 6013 | BE 2030 | ■ |
| 10mm | BE 6010 | BE 2018 | ■ | 25mm | BE 6014 | BE 2033 | ■ |
| 12mm | LP 16-012 | LP 20-012 | ■ | 28mm | BE 6015 | — | ■ |
| 14mm | BE 6011 | BE 2022 | ■ | | | | |

Blom-Singer® Special Order Low Pressure Voice Prostheses

The Blom-Singer® Special Order Low Pressure is a patient changeable, one-way valve prosthesis.


| Length | 16 Fr Special Increased Resistance | 16 Fr Special Length | 20 Fr Special Length | 20 Fr Special Increased Resistance | Rx |
|--------|------------------------------------|----------------------|----------------------|------------------------------------|----|
| 4mm | — | SPECIAL 11-902-06 | SPECIAL 11-251-06 | — | ■ |
| 6mm | LP 16-006-IR | — | — | LP 20-006-IR | ■ |
| 8mm | LP 16-008-IR | — | — | LP 20-008-IR | ■ |
| 10mm | LP 16-010-IR | — | — | LP 20-010-IR | ■ |
| 12mm | LP 16-012-IR | — | — | LP 20-012-IR | ■ |
| 14mm | LP 16-014-IR | — | — | LP 20-014-IR | ■ |
| 16mm | — | SPECIAL 11-902-04 | SPECIAL 11-251-04 | — | ■ |
| 18mm | LP 16-018-IR | — | — | — | ■ |
| 20mm | — | SPECIAL 11-902-05 | SPECIAL 11-251-05 | — | ■ |
| 22mm | LP 16-022-IR | — | — | — | ■ |
| 25mm | LP 16-025-IR | — | — | — | ■ |
| 28mm | LP 16-028-IR | — | — | — | ■ |

Blom-Singer® Voice Prosthesis Sizers

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

The Blom-Singer Voice Prosthesis Sizer is a device used by clinicians to measure the length of the tracheoesophageal puncture to assist in selection of the appropriate voice prosthesis.

| Ref | Description |
|---------|---|
| BE 6036 | Voice Prosthesis Sizer with 3 inserters (for 16 Fr voice prostheses) |
| BE 2036 | Voice Prosthesis Sizer with 3 inserters (for 20 Fr voice prostheses) |
| BE 6096 | Replacement Inserters (for use with both BE 6036 and BE 2036 Sizers), 10 ea |

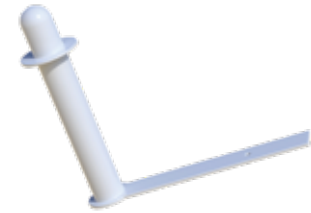

Blom-Singer® Dilator-Sizers

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

The Blom-Singer Dilator-Sizer combines two tasks into one innovative product: sizing and dilation. In one step, the clinician can dilate and measure the tracheoesophageal puncture during voice prosthesis placement procedures.

| Ref | Description |
|---------|--|
| BE 2051 | 18 Fr Dilator-Sizer (for 16 Fr prostheses) |
| BE 2052 | 22 Fr Dilator-Sizer (for 20 Fr prostheses) |

Blom-Singer® Tracheoesophageal Puncture Dilators

| Ref | Description | Rx |
|---------|---|----|
| BE 6050 | 18 Fr Puncture Dilator-Sizer (for 16 Fr prostheses) | ■ |
| BE 2050 | 22 Fr Puncture Dilator-Sizer (for 20 Fr prostheses) | ■ |

VOICE PROSTHESIS ACCESSORIES

Blom-Singer® Replacement Gel Caps

Use of the Blom-Singer Replacement Gel Caps provides a smooth, rounded shape to the esophageal flange end of the voice prosthesis, enabling easy insertion into the tracheoesophageal puncture.

| Ref | Description | Rx |
|---------|------------------------------------|----|
| BE 3190 | Replacement Gel Caps, 16 Fr, 90 ea | ■ |
| BE 3195 | Replacement Gel Caps, 18 Fr, 90 ea | ■ |
| BE 3290 | Replacement Gel Caps, 20 Fr, 90 ea | ■ |
| BE 3295 | Replacement Gel Caps, 22 Fr, 60 ea | ■ |

Blom-Singer® Flushing Device and Cleaning Brushes

For routine cleaning of a voice prosthesis as recommended by your clinician.

| Ref | Description |
|----------|--|
| IN 4050 | Flushing Device, pack of 3 |
| IN 4064* | Cleaning Brush, 16/20 Fr, Small, 4mm–8mm, 3 ea |
| IN 4065* | Cleaning Brush, 16/20 Fr, Medium, 9mm–16mm, 3 ea |
| IN 4066* | Cleaning Brush, 16/20 Fr, Large, 18mm–22mm, 3 ea |

*Fits indwelling style voice prostheses only.

Blom-Singer® Plug Inserts

The Blom-Singer Plug Insert is designed for temporary management of leakage through a voice prosthesis. Not compatible with patient changeable voice prostheses.

| Ref | Description | Rx |
|---------|--|----|
| IN 4069 | 16 Fr Plug Insert (Advantage and Classic, 6mm to 20mm) | ■ |
| IN 4070 | 20 Fr Plug Insert (Advantage and Classic, 6mm to 20mm) | ■ |



Blom-Singer® HMEs

Simple choices to fit your lifestyle

Left: Joanne
Year of TL: 2015
Wearing AccuFit
housing with
SpeakFree HME
Hands-Free

Right: Meagon
Year of TL: 2016
Wearing HydroFit
housing with
Day&Night HME

Day&Night® HME



- Wear all day and all night, including bedtime
- Attaches to the entire family of InHealth baseplates and StomaSoft tubes.



Day&Night®
EasyFlow® HME and
ClassicFlow® HME

SpeakFree® HME

- Your HME and hands free speaking valve, all in one.
- Attaches to the entire family of InHealth baseplates and StomaSoft tubes.



SpeakFree®
EasyFlow® HME and
ClassicFlow® HME



HEAT & MOISTURE EXCHANGE (HME)

Research shows that long-term HME use results in less mucus¹, less coughing², and decreased shortness of breath, fatigue, sleep problems, anxiety, and depression^{1,2}.

Blom-Singer® SpeakFree® Hands-Free Valve with ClassicFlow® HME

Blom-Singer SpeakFree HME Hands-Free Valve is the first heat and moisture exchange (HME) cartridge with a single-use fully integrated adjustable hands-free valve. The SpeakFree HME is translucent and non-reflective for a natural look. The SpeakFree ClassicFlow HME is designed for increased moisture retention and resistance.

| Ref | Description |
|-----------|--|
| BE 1090EZ | SpeakFree Hands-Free Valve with ClassicFlow HME, 30 per pack |



SpeakFree
ClassicFlow

Blom-Singer® SpeakFree® Hands-Free Valve with EasyFlow® HME

Blom-Singer SpeakFree HME Hands-Free Valve is the first heat and moisture exchange (HME) cartridge with a single-use fully integrated adjustable hands-free valve. The SpeakFree HME is translucent and non-reflective for a natural look. With SpeakFree EasyFlow, the foam filter is more porous and allows you to breathe more freely to match your activity level and pulmonary needs.

| Ref | Description |
|-----------|---|
| BE 1090EF | SpeakFree Hands-Free Valve with EasyFlow HME, 30 per pack |



SpeakFree
EasyFlow

Blom-Singer® Day&Night® ClassicFlow® HME

The Blom-Singer Day&Night ClassicFlow can be worn all day and night including bedtime. The Day&Night HME is translucent and non-reflective for a natural look.

| Ref | Description |
|---------|---|
| BE 1091 | Day&Night ClassicFlow HME Cartridge, pack of 30 (24-hour use) |



Day&Night
ClassicFlow

Blom-Singer® Day&Night® EasyFlow® HME

The Blom-Singer Day&Night EasyFlow allows you to breathe more freely to match your activity level and pulmonary needs. The Day&Night HME is translucent and non-reflective for a natural look.

| Ref | Description |
|---------|--|
| BE 1092 | Day&Night EasyFlow HME Cartridge, pack of 30 (24-hour use) |



Day&Night
EasyFlow

1. Ackerstaff AH, Hilgers FJ, et al. Improvements in respiratory and psychosocial functioning following total laryngectomy by the use of a heat and moisture exchanger. *Annals of Otolaryngology, Rhinology & Laryngology*. 1993;102(11):878-883.

2. Ackerstaff AH, Hilgers FJ, et al. Heat and moisture exchangers as a treatment option in the post-operative rehabilitation of laryngectomized patients. *Clinical Otolaryngology*. 1995;20:504-509.

ADHESIVE HOUSINGS

InHealth Technologies® presents the FitSeries—adhesive housings designed to meet your skin and stoma needs after a total laryngectomy. Available in two styles, three shapes each. Compatible with Day&Night and SpeakFree Hands-Free HME cartridges.

Blom-Singer® AccuFit® Adhesive Housings


The Blom-Singer AccuFit Adhesive Housing is designed to be our strongest adhesive option while still maintaining flexibility. It has increased elasticity and is less prone to tearing. Available in three shapes: round, oval, and oval extra.

| Ref | Description |
|---------|--|
| BE 6082 | AccuFit Adhesive Housing, Round, pack of 30 |
| BE 6083 | AccuFit Adhesive Housing, Oval, pack of 30 |
| BE 6084 | AccuFit Adhesive Housing, Oval Extra, pack of 30 |

Blom-Singer® HydroFit® Adhesive Housings

The Blom-Singer HydroFit Adhesive Housing is a hydrocolloid adhesive baseplate made with a protective gel-like layer and a moisture-absorbent adhesive. The HydroFit is gentle on the skin and flexible, but maintains strong adhesion. Available in three shapes: round, oval, and oval extra.

| Ref | Description |
|---------|---|
| BE 6085 | HydroFit Adhesive Housing, Round, pack of 30 |
| BE 6086 | HydroFit Adhesive Housing, Oval, pack of 30 |
| BE 6087 | HydroFit Adhesive Housing, Oval Extra, pack of 30 |

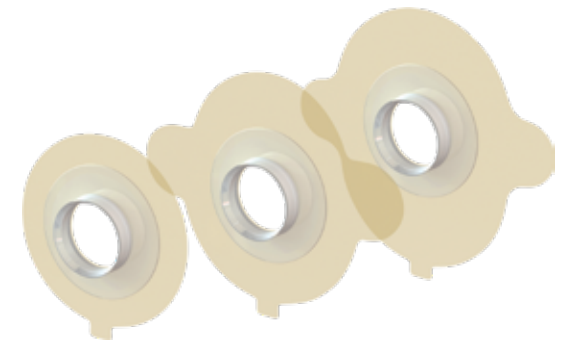

TRACHEOSTOMA ATTACHMENTS

Blom-Singer® Tracheostoma Valve Housings

The Tracheostoma Valve Housing is a reusable, flexible housing that attaches to the Day&Night and SpeakFree Hands-Free HME cartridge or Shower Guard over the stoma. It is applied with double sided tape discs. It is available in two sizes, standard and large.

| Ref | Description |
|---------|--|
| BE 6038 | Tracheostoma Valve Housing, standard, 1 ea |
| BE 6039 | Tracheostoma Valve Housing, large, 1 ea |




Blom-Singer® Adhesive Tape Discs

A variety of double-sided tape discs are available to provide attachment of the tracheostoma valve housing.

| Ref | Description | Ref | Description |
|---------|--|---------|---|
| BE 6041 | Regular Tape Discs, standard, pack of 30 | BE 6044 | Thin Tape Discs, large, pack of 30 |
| BE 6042 | Regular Tape Discs, large, pack of 30 | BE 6034 | Heavy Duty Tape Discs, standard, pack of 30 |
| BE 6043 | Thin Tape Discs, standard, pack of 30 | BE 6035 | Heavy Duty Tape Discs, large, pack of 30 |



Kevin
Year of TL: 2001
Wearing SpeakFree
Hands-Free HME

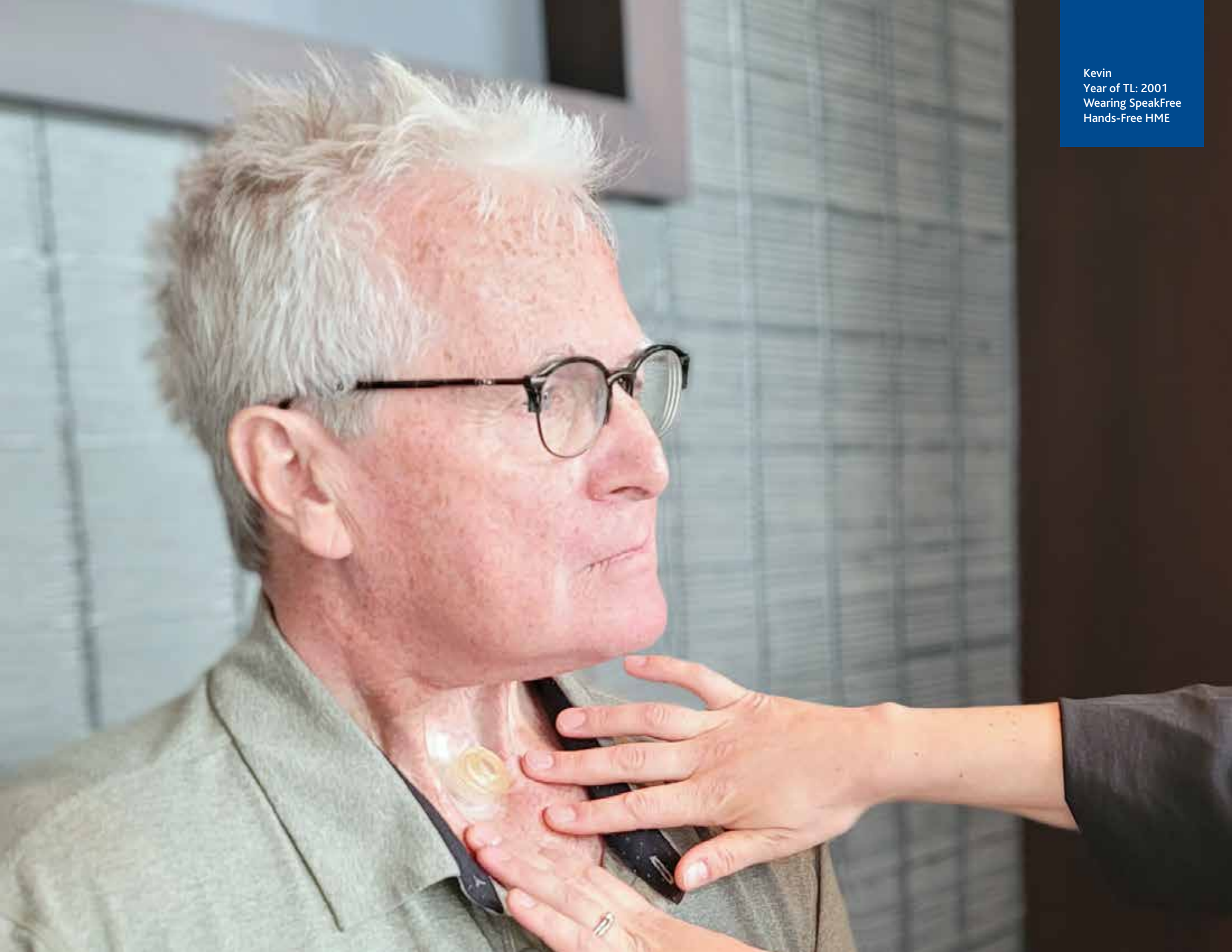

● AND BEYOND

HANDS-FREE SPEECH VALVE

Blom-Singer® SpeakFree® Hands-Free Valve with ClassicFlow® HME

Blom-Singer SpeakFree HME Hands-Free Valve is the first heat and moisture exchange (HME) cartridge with a single-use fully integrated adjustable hands-free valve. The SpeakFree HME is translucent and non-reflective for a natural look. The SpeakFree ClassicFlow HME is designed for increased moisture retention and resistance.

| Ref | Description |
|-----------|--|
| BE 1090EZ | SpeakFree Hands-Free Valve with ClassicFlow HME, 30 per pack |

Blom-Singer® SpeakFree® Hands-Free Valve with EasyFlow® HME

Blom-Singer SpeakFree HME Hands-Free Valve is the first heat and moisture exchange (HME) cartridge with a single-use fully integrated adjustable hands-free valve. The SpeakFree HME is translucent and non-reflective for a natural look. With SpeakFree EasyFlow, the foam filter is more porous and allows you to breathe more freely to match your activity level and pulmonary needs.

| Ref | Description |
|-----------|---|
| BE 1090EF | SpeakFree Hands-Free Valve with EasyFlow HME, 30 per pack |



SpeakFree
ClassicFlow



SpeakFree
EasyFlow

SKIN CARE

Skin-Prep™ Liquid Film-Forming Protective Barrier Wipes

Skin-Prep is a liquid film-forming dressing that, upon application to intact skin, forms a protective film to help reduce friction during removal of tapes and films.

| Ref | Description |
|--------|---|
| 420400 | Skin-Prep Liquid Film-Forming Protective Barrier Wipes, box of 50 |

Skin Tac™ Adhesive Barrier Wipes

Skin Tac liquid adhesive barrier wipes are latex-free, hypo-allergenic, clear, non-rubber adhesive barrier wipe.

| Ref | Description |
|---------|-----------------------------------|
| MS407W* | Adhesive Barrier Wipes, box of 50 |

*May not be available in all international markets.

NFS Adhesive (Silicone Adhesive)

PSA Pressure Sensitive Adhesive, silicone based adhesive for use with all silicone prosthetic devices.

| Ref | Description |
|-------|--------------------|
| B-530 | NFS Adhesive, 1 oz |

Remove™ Adhesive Remover

Gently cleans hydrocolloid, acrylic and rubber-based adhesive residues from the skin without irritation, discomfort or trauma.

| Ref | Description |
|----------|-----------------------------|
| 59403125 | Adhesive Remover, box of 50 |



STOMA PROTECTION AND ACCESSORIES

ADDvox® Stoma Cover Mock Turtlenecks (Only available in the US)

The ADDvox Stoma Cover Mock Turtleneck has inner layers of material which helps keep debris from entering the stoma. They slip under a shirt like a fashion dickey and have an outer layer of soft 100% cotton knit. Each cover is 9-1/2" long x 8-3/4" wide. Velcro fastener adjusts to fit necks of 12-1/2" to 22". Each color sold separately.

| Ref | Description | Color | Quantity |
|------------|-----------------------------|-------|----------|
| BE 6227-R2 | Stoma Cover Mock Turtleneck | White | 1 ea |
| BE 6228-R2 | Stoma Cover Mock Turtleneck | Black | 1 ea |
| BE 6229-R2 | Stoma Cover Mock Turtleneck | Navy | 1 ea |



ADDvox® Stoma Scarves (Only available in the US)

The ADDvox Stoma Cover Scarf has inside layers to help filter out dust, dirt and pollen, is completely washable, and needs little care. This attractive neckline scarf looks just like a regular scarf or ascot and is available in a wide range of colors for men or women. All scarves have Velcro closures and adjust to fit necks from 17-3/4" to 22" (can be made smaller with tucks or stitching). Each color sold separately.

| Ref | Description | Color | Quantity |
|------------|-------------------|-------|----------|
| BE 6232-R2 | Stoma Cover Scarf | Black | 1 ea |
| BE 6233-R2 | Stoma Cover Scarf | White | 1 ea |
| BE 6234-R2 | Stoma Cover Scarf | Navy | 1 ea |



ADDvox® Cloth Bib-Style Stoma Covers (Only available in the US)

The ADDvox soft, long-lasting bib protects the tracheostoma and lets you breathe more easily. This "bib-style" Stoma Cover will look crisp and clean after repeated laundering. Five soft layers keep dust and dirt from reaching the trachea. Fits necks up to 22". Velcro closure. Each color sold separately.

| Ref | Description | Color | Quantity |
|------------|-----------------------------|-------|----------|
| BE 6224-R2 | Cloth Bib-Style Stoma Cover | White | 1 ea |
| BE 6225-R2 | Cloth Bib-Style Stoma Cover | Beige | 1 ea |
| BE 6226-R2 | Cloth Bib-Style Stoma Cover | Navy | 1 ea |



ADDvox® Silicone Shower Shield (Only available in the US)

The one-piece ADDvox Silicone Shower Shield is worn during showering to deflect the water away from your stoma. Extremely soft and pliable, it creates a soft, comfortable and conforming seal along the irregular contours of your neck, and vents on the underside allowing for airflow. Made of medical grade silicone, this latex-free Silicone Shower Shield is comfortable, durable, and not affected by temperature, ozone, sunlight, soap or detergents.

| Ref | Description |
|------------|------------------------------|
| BE 6079-R2 | Silicone Shower Shield, 1 ea |

Blom-Singer® Shower Guard

The Blom-Singer Shower Guard is a one-piece, right-angle plastic shower guard that comes with a Tracheostoma Standard Adhesive Housing and 30 double-sided Adhesive Tape Discs. Also compatible with the Blom-Singer StomaSoft® Laryngectomy Tube.

| Ref | Description |
|---------|--|
| BE 6048 | Shower Guard—Tracheostoma Valve Housing (BE 6038) and Regular Tape Discs, pk 30 (BE 6041) included |



Blom-Singer® Foam Stoma Protector

A simple and temporary (disposable) solution to help prevent debris from entering the stoma. Adhesive strip on one edge allows for easy attachment, access to the stoma to occlude for speech with a voice prosthesis, or for coughing. Measurement listed is approximate.

| Ref | Description |
|---------|--|
| BE 6220 | Foam Stoma Protector, Tan, Large 2.7" L x 2.7" W (7cm x 7cm), pack of 30 |



ADDvox® Stoma Filters with Microporous Adhesive

ADDvox Stoma Filters offer protection against dust, dirt and pollen. Designed for single use with a hypoallergenic, non-irritating, self-adhesive microporous strip that is comfortable and easy to apply. Not recommended for cannula users. 2-1/4" wide x 2-1/2" high x 1/8" thick. Available in off-white and skin tone.

| Ref | Description |
|------------|--------------------------------------|
| BE 6222-R2 | Stoma Filters, Skin tone, pack of 60 |
| BE 6223-R2 | Stoma Filters, Off-white, pack of 60 |





Professional Products

VOICE PROSTHESIS MANAGEMENT

Blom-Singer® Voice Prosthesis Sizers

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

The Blom-Singer Voice Prosthesis Sizer is a device used by clinicians to measure the length of the tracheoesophageal puncture to assist in selection of the appropriate voice prosthesis.

| Ref | Description |
|---------|---|
| BE 6036 | Voice Prosthesis Sizer with 3 inserters (for 16 Fr voice prostheses) |
| BE 2036 | Voice Prosthesis Sizer with 3 inserters (for 20 Fr voice prostheses) |
| BE 6096 | Replacement Inserters (for use with both BE 6036 and BE 2036 Sizers), 10 ea |



Blom-Singer® Dilator-Sizers

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

The Blom-Singer Dilator-Sizer combines two tasks into one innovative product: sizing and dilation. In one step, the clinician can dilate and measure the tracheoesophageal puncture during voice prosthesis placement procedures.

| Ref | Description |
|---------|--|
| BE 2051 | 18 Fr Dilator-Sizer (for 16 Fr prostheses) |
| BE 2052 | 22 Fr Dilator-Sizer (for 20 Fr prostheses) |



Blom-Singer® Replacement Gel Caps

| Ref | Description | Rx |
|---------|------------------------------------|----|
| BE 3190 | Replacement Gel Caps, 16 Fr, 90 ea | ■ |
| BE 3195 | Replacement Gel Caps, 18 Fr, 90 ea | ■ |
| BE 3290 | Replacement Gel Caps, 20 Fr, 90 ea | ■ |
| BE 3295 | Replacement Gel Caps, 22 Fr, 60 ea | ■ |



Manometer with Manometer Adaptor

A tool to teach patients how to self-regulate intratracheal air pressure. Proper airflow control may allow for better voice quality, less fatigue, and longer duration of the tracheostoma adhesive housing seal.

| Ref | Description |
|---------|---|
| BE 8000 | Manometer Kit includes 1 manometer, 1 manometer adaptor and 1 adhesive housing |
| BE 8020 | Replacement Manometer Adaptor includes 5 manometer adaptors and 5 adhesive housings |




Blom-Singer® Insufflation Test Set

The Blom-Singer Insufflation Test Set is designed to aid the clinician in subjectively predicting the ability of the laryngectomized patient to produce tracheoesophageal voice.

| Ref | Description |
|---------|---------------------------------|
| BE 6040 | Insufflation Test Set, box of 6 |



Guiermo
Year of TL: 2021
Wearing Accufit housing
and Day&Night HME

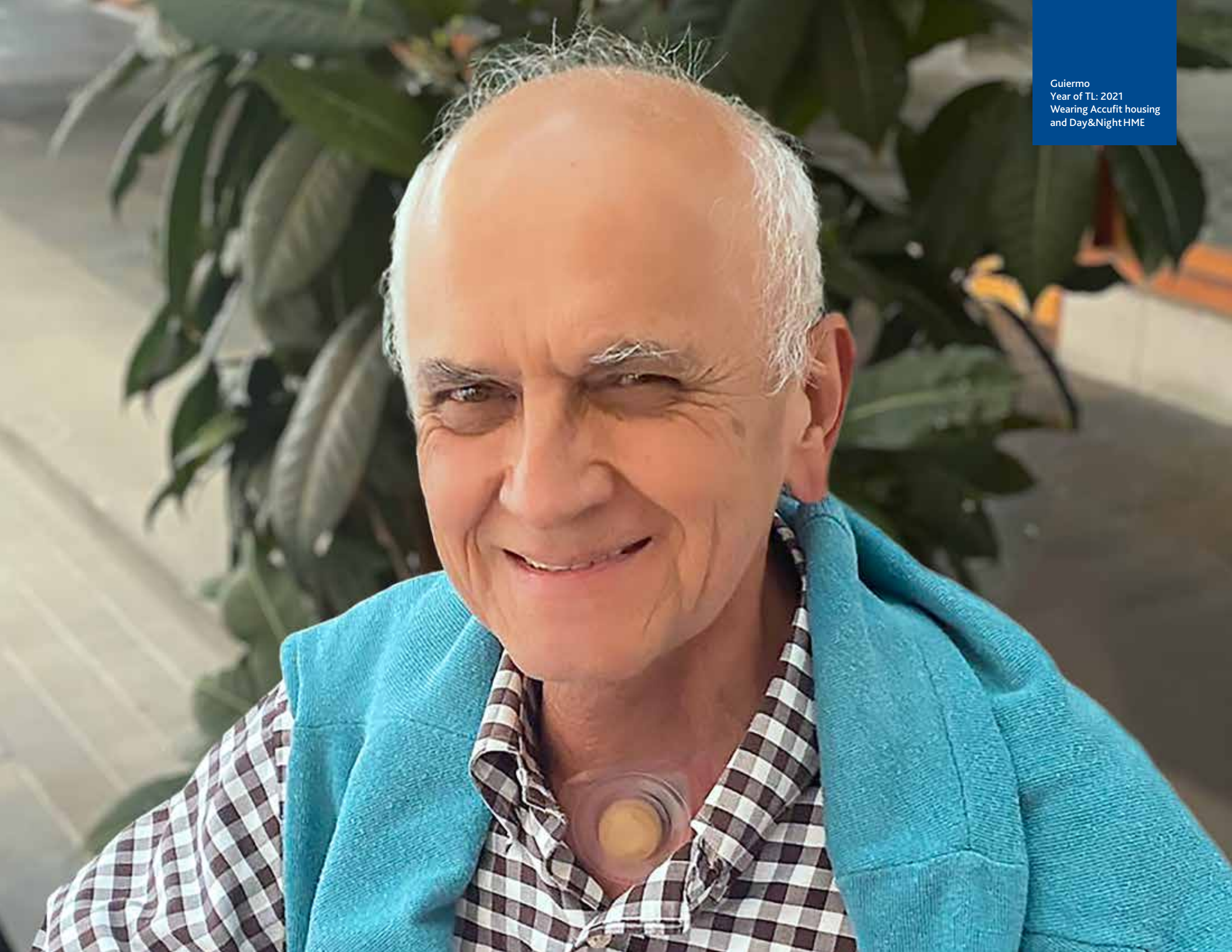

Katz Extractor[®]

OTO-RHINO FOREIGN BODY REMOVER

The easy way out



Katz Extractor® Oto-Rhino Foreign Body Remover

Designed in conjunction with David L. Katz, MD, MPH, the Katz Extractor Oto-Rhino Foreign Body Remover is intended for extraction of impacted foreign bodies in the nasal passage. Single use only. This device may also be used in the external ear canal, based upon clinical judgment. Covered by one or more of the following patents and other patents pending: US#5,454,817. Katz Extractor® is a registered trademark of Freudenberg Medical, LLC.

| Ref | Description |
|---------|--|
| FB 5000 | Oto-Rhino Foreign Body Remover, box of 5 |

TO LEARN MORE
Scan QR code
to visit the
Katz Extractor website





Nasal Septal Perforation Prosthesis

A two-piece, magnetically coupled solution for non-surgical closure of nasal septal perforations.

Nasal Septal Perforation Prostheses

The Blom-Singer® Nasal Septal Perforation Prosthesis is a two-piece, magnetically coupled solution for non-surgical closure of nasal septal perforations.

OVAL

| Ref | Description |
|---------|--|
| SP 2112 | Nasal Septal Perforation Prosthesis Oval 11 x 13 x 2mm |
| SP 2113 | Nasal Septal Perforation Prosthesis Oval 11 x 13 x 3mm |
| SP 2132 | Nasal Septal Perforation Prosthesis Oval 13 x 17 x 2mm |
| SP 2133 | Nasal Septal Perforation Prosthesis Oval 13 x 17 x 3mm |
| SP 2152 | Nasal Septal Perforation Prosthesis Oval 15 x 21 x 2mm |
| SP 2153 | Nasal Septal Perforation Prosthesis Oval 15 x 21 x 3mm |









Diagrams are NOT actual size.
Flange width is 5mm for all prostheses.

ROUND

| Ref | Description |
|---------|--|
| SP 1052 | Nasal Septal Perforation Prosthesis Round 5 x 2mm |
| SP 1053 | Nasal Septal Perforation Prosthesis Round 5 x 3mm |
| SP 1072 | Nasal Septal Perforation Prosthesis Round 7 x 2mm |
| SP 1073 | Nasal Septal Perforation Prosthesis Round 7 x 3mm |
| SP 1102 | Nasal Septal Perforation Prosthesis Round 10 x 2mm |
| SP 1103 | Nasal Septal Perforation Prosthesis Round 10 x 3mm |



Diagrams are NOT actual size.
Flange width is 5mm for all prostheses.

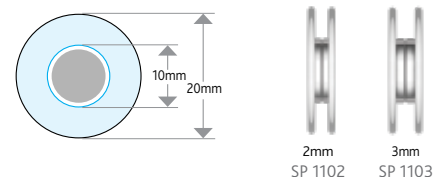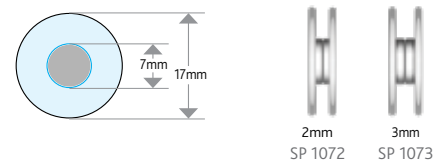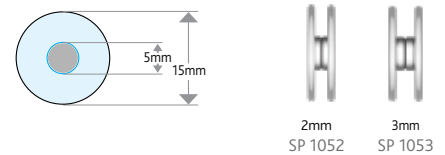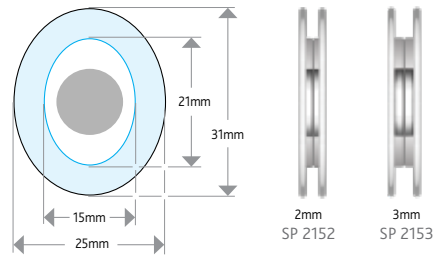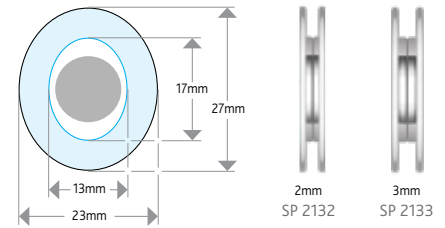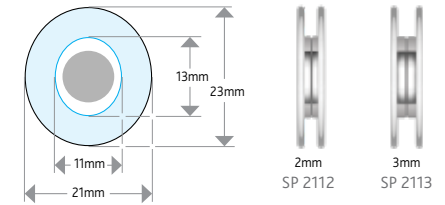

TO LEARN MORE
Scan QR code to visit the
Nasal Septal Perforation Prostheses website

ORDERING INFORMATION

UNITED STATES

Phone: Toll-Free: 800.477.5969

Fax: Toll-Free: 888.371.1530

E-mail: order@inhealth.com

Prescription E-mail: ihtrx@inhealth.com

Order Online: www.inhealth.com

Mail:

InHealth Technologies®

1110 Mark Ave.

Carpinteria, CA 93013 • USA

CANADA

Phone: Toll-Free: 1.888.437.4123

Fax: Toll-Free: 1.226.641.5285

E-mail: canadaorder@inhealth.com

Website: www.inhealth.com/canada

Mail:

InHealth Canada

65 Spruce Street

Tillsonburg ON N4G 5C4

**Rx: Products
with the Rx
notation require
a prescription
from a physician
(US patients).**

NEW PATIENT ENROLLMENT INSTRUCTIONS

Enrolling as a new patient? It's easy! Here is what to do:

1. Obtain a **prescription** from your provider.
2. Fill out the **Patient Enrollment Form**. All forms are available under the **Resources** tab on the InHealth Technologies® website: www.inhealth.com
3. Send back to InHealth Technologies® with a copy of your insurance card (front and back). You can send via email: insurance@inhealth.com or fax: **888.371.1530**
4. Place your order! Please remember that a prescription alone will not automatically generate an order.

Need insurance or Medicare help? Contact us!

Phone: **800.477.5969 option 3**

Email: medicare@inhealth.com

While you don't have to fill anything out on the **Patient Bill of Rights** form, please make sure you review it. We want you to know your rights.

FREE EDUCATIONAL RESOURCES

- Available on www.inhealth.com under the **Resources** tab.
- Visit InHealth Technologies® video channel at <https://inhealth.com/videos> to view our videos

In 1978, Drs. Singer and Blom pioneered the methods of tracheoesophageal puncture and valved silicone voice prostheses that have been the international standard for voice restoration for over 45 years, allowing thousands of patients worldwide to regain their ability to speak. As a leading supplier of advanced voice restoration systems, our goal at InHealth Technologies® is to help laryngectomees regain and preserve a vital part of their lives—their voice. We are proud to continue the Blom-Singer® legacy today with clinical support, input, and leadership.

Blom-Singer®





Janice
Year of TL: 2014
Wearing Accufit housing
with SpeakFree HME
Hands-Free Valve

STAY CONNECTED...



PATIENTS
SCAN QR CODE

Speakers Club™

Patients: Scan the QR code to sign up today!

We created this digital E-News for you! **Speakers Club E-News** will provide you with information on upcoming patient events, new products and resources.



PROFESSIONALS
SCAN QR CODE

ProNET™

Professionals: Scan the QR code to sign up today!

We created this digital E-News specifically for professionals. **ProNET E-News** will deliver useful laryngectomy and voice prosthesis information for you and your practice, including material you may want to share with your patients.

INHEALTH
TECHNOLOGIES®

United States
1110 Mark Ave., Carpinteria, CA 93013
Toll-Free: 800.477.5969
Toll-Free Fax: 888.371.1530

International
Tel: +1.805.684.9337
Fax: +1.805.684.8594
E-mail: international@inhealth.com

Terms & Conditions: <https://inhealth.com/terms-conditions>
©2025 InHealth Technologies®—Manufactured by Freudenberg Medical, LLC (241104.00)

Blom-Singer®

The voice of experience since 1978



Order online at www.inhealth.com