


INTRODUCING  
**FirstFit™ Surgical Kit**  
with Voice Prosthesis ready for insertion

# A Better Solution for TE Puncture\*


## FirstFit™ Surgical Kit

The FirstFit™ Surgical Kit is designed to create a primary or secondary tracheoesophageal puncture, to facilitate the placement of the included sterile Blom-Singer® Classic™ Voice Prosthesis.

Everything you need, in one sterile package. Ergonomic design of the needle and pharynx protector allowing for more control and accuracy. Color contrast of the guidewire and inserter allowing for easy detection to adjust grip and/or slow down for steady delivery of the voice prosthesis.

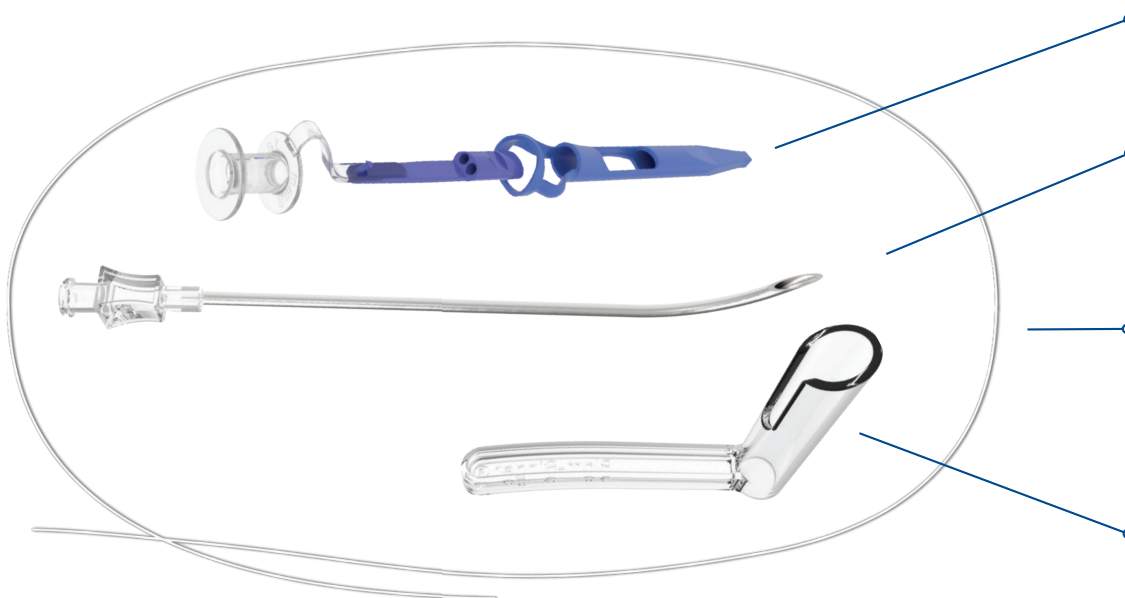

- **Inserter:**  
Retrograde insertion with pre-loaded voice prosthesis
- **13 Gauge Surgical Stainless Steel Needle:**  
New ergonomic curve of needle to allow for more control and accuracy at puncture site
- **91cm Guidewire:**  
Reversible for easy correction. Color contrast of the guidewire and inserter allowing for easy detection to adjust grip and/or slow down for steady delivery
- **Pharynx Protector:**  
The shape, size and transparency is more universally user-friendly

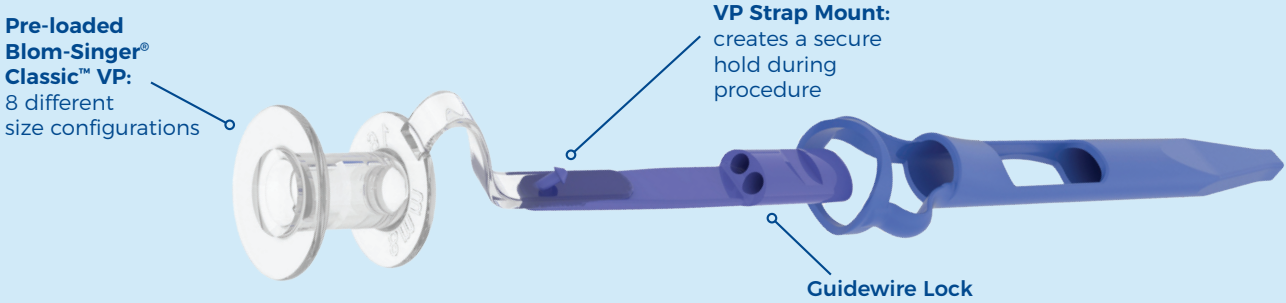

One sterile package, packed for ease of use



Retrograde  
insertion



One sterile package  
with pre-loaded  
Blom-Singer® voice  
prosthesis



Color contrast in tools  
to allow more visual  
detection



Ergonomic design  
for more control  
and accuracy

For more product information, visit [www.firstfitsurgicalkit.com](http://www.firstfitsurgicalkit.com)

OPTIONS TO FIT ALL OF YOUR PATIENTS

Ref	Description
TP 1608	Sterile FirstFit Surgical Kit + Classic 16 Fr 8 mm Voice Prosthesis*
TP 1610	Sterile FirstFit Surgical Kit + Classic 16 Fr 10 mm Voice Prosthesis*
TP 1612	Sterile FirstFit Surgical Kit + Classic 16 Fr 12 mm Voice Prosthesis*
TP 1614	Sterile FirstFit Surgical Kit + Classic 16 Fr 14 mm Voice Prosthesis*
TP 2008	Sterile FirstFit Surgical Kit + Classic 20 Fr 8 mm Voice Prosthesis*
TP 2010	Sterile FirstFit Surgical Kit + Classic 20 Fr 10 mm Voice Prosthesis*
TP 2012	Sterile FirstFit Surgical Kit + Classic 20 Fr 12 mm Voice Prosthesis*
TP 2014	Sterile FirstFit Surgical Kit + Classic 20 Fr 14 mm Voice Prosthesis*

\*This product is restricted to sale or use by a licensed professional only.