



# Katz Extractor<sup>®</sup>

OTO-RHINO FOREIGN BODY REMOVER

## The easy way out

The fast, flexible way to remove foreign bodies from noses and ears

- **Efficient and Safe**  
Great for restless patients.
- **Reliable**  
Trusted by healthcare professionals worldwide.
- **Convenient**  
Lightweight, fully assembled and ready to use.

**Order the  
Katz Extractor<sup>®</sup>  
Five-pack today!**

# Katz Extractor®

OTO-RHINO FOREIGN BODY REMOVER

A quick and simple solution for foreign body removal from the nose or ear.




The Katz Extractor® Oto-Rhino Foreign Body Remover is a small, lightweight, single-use device for quick and easy extraction of foreign bodies from the nose or external ear canal. It comes ready to use in a five-pack dispenser box. There are no expensive supplies, such as forceps, hooks or balloon tip catheters/syringes to locate and assemble. No additional suction equipment is required.

To see a demonstration video, visit: [www.katzextractor.com](http://www.katzextractor.com)




Please Scan  
to access website  
and IFU

*It takes just three simple steps to complete an extraction with the Katz Extractor!*

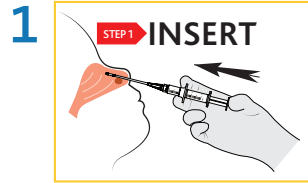

Insert past foreign body

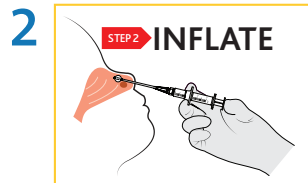

Depress the syringe plunger to inflate balloon



Extract catheter out of canal with balloon inflated

ORDER INFORMATION: Product Number: FB5000

To order online visit:  
[www.katzextractor.com](http://www.katzextractor.com)

Tel: 800.477.5969

Email: [order@inhealth.com](mailto:order@inhealth.com)

#### CUSTOMER SERVICE

Our knowledgeable representatives are available to assist you  
Monday through Thursday 6:30 am to 4:00 pm  
and Friday 6:30 am to 2:00 pm, PT.

#### SHIPPING

All orders ship via our standard five-business-day service.  
Overnight shipping is available upon request for orders placed before 2:00 pm, PT.

#### Device Coding and Reimbursement

CPT Code	Procedure Description	Reimbursement		
		Ambulatory Surgical Center	Outpatient Surgery at Hospital	Physician Fee Schedule
30300	Removal foreign body, intranasal	Facility Payment (National Medicare Average) \$570.23	APC (Ambulatory Payment Classifications) \$107.10	Payment (National Medicare Average) \$201.20
69200	Removal foreign body from external auditory canal; w/o general anesthesia	\$522.35	\$107.10	\$56.02

Designed in conjunction with David L. Katz, M.D., MPH, the Katz Extractor is intended principally for extraction of impacted foreign bodies in the nasal passages. This device may also be used in the external ear canal, based upon clinical judgement.

**INHEALTH**  
TECHNOLOGIES®

Stay Connected. Like, Share, Follow.



©2025 InHealth Technologies

Manufactured by Freudenberg Medical, LLC (240806.02)