


**INHEALTH**  
TECHNOLOGIES®

2024 PRODUCT CATALOG

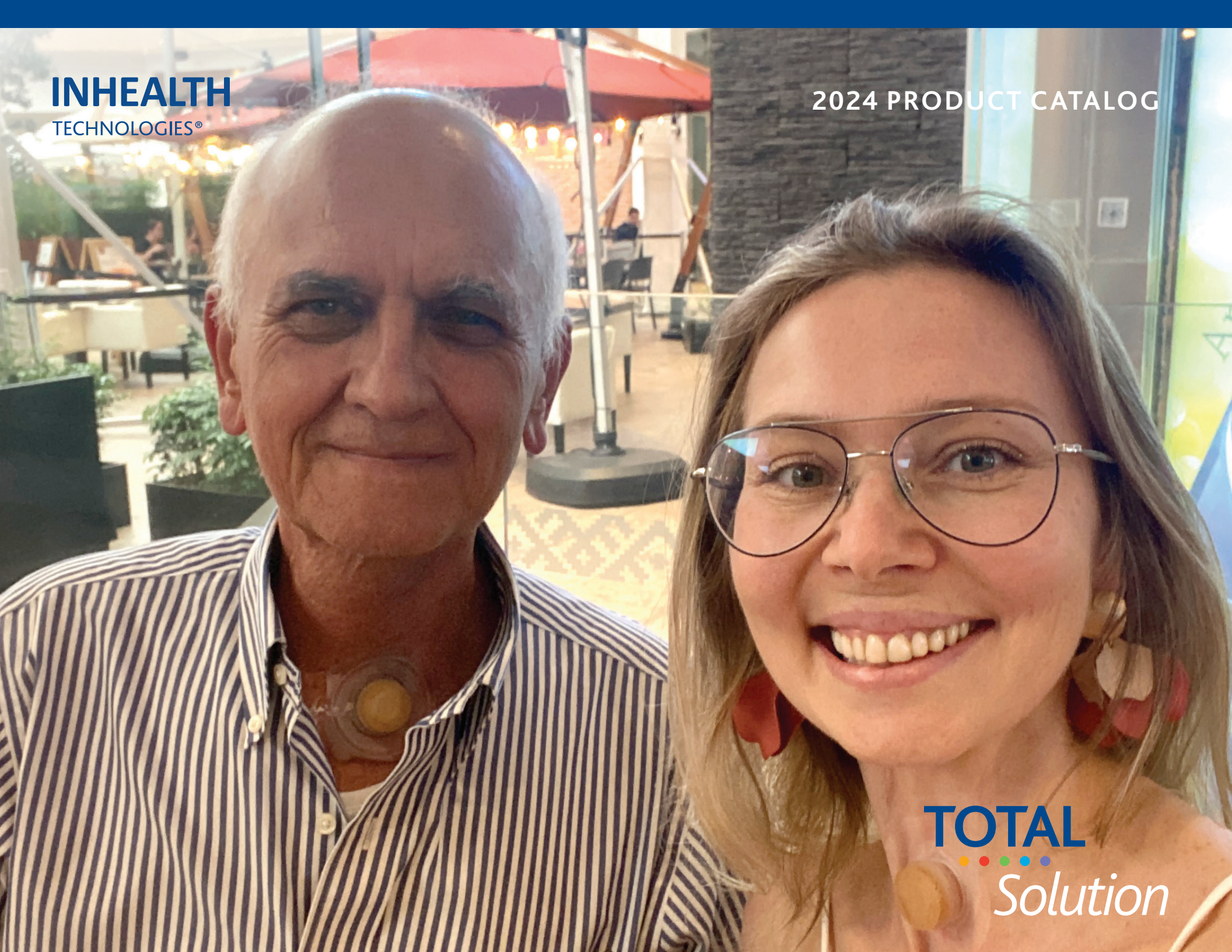

**TOTAL**  
  
*Solution*

# CONTENTS

<b>Pre-Operative Education</b>	<b>2</b>
"What to Expect" Brochure	2
Educational Videos	2
<b>Surgical Solutions</b>	<b>3</b>
Tracheoesophageal Puncture Set (Sterile)	3
Indwelling Voice Prosthesis (Sterile)	3
Adjustable Bi-Flanged Fistula Prosthesis (Sterile)	3
<b>Post-Operative Care</b>	<b>4</b>
Laryngectomy Tubes	4
Marpac Tracheostomy Tube Holder/Collar	4
StomaSoft Cleaning Swabs	4
Electrolarynx	5
Post-Operative Care Kits (For Total Laryngectomy)	6
<b>Next Steps</b>	<b>8</b>
Voice Prosthesis (Clinician Placed)	8
Special Order Voice Prosthesis (Clinician Placed)	10
Voice Prosthesis (Patient Changeable)	12
Voice Prosthesis Accessories	13
Heat & Moisture Exchange (HME)	14
Adhesive Housings	16
Tracheostoma Attachments	16
<b>And Beyond</b>	<b>18</b>
Hands-Free Speech Valves	18
Skin Care	18
Stoma Protection and Accessories	19
<b>Professional Products</b>	<b>21</b>
Voice Prosthesis Management	22
RENÚ® Voice and RENÚ® Gel	23
Oto-Rhino Foreign Body Remover	24
Nasal Septal Perforation Prosthesis	26
<b>Ordering Information</b>	<b>27</b>
<b>New Patient Enrollment Instructions</b>	<b>27</b>
<b>Free Educational Resources</b>	<b>27</b>

## InHealth Technologies® offers a TOTAL Solution for your total laryngectomy from pre- to post-op and beyond

TOTAL Solution serves patients at all phases of voice restoration, from products that deliver solutions for pre- and post-operative care and beyond. This continuum of care is enhanced by InHealth Technologies' customer-focused team, including leading clinical specialists, experienced patient consultants and dedicated customer service members.

Cover Image:

Left: Guiermo  
Year of TL: 2021  
Wearing Accufit housing  
and Day&Night HME

Right: Ksenia  
Year of TL: 2019  
Wearing Accufit housing  
and Day&Night HME

# TOTAL Solution

## Pre-Operative Education

- "What to Expect" Brochure
- Patient Videos

PHASE

1

PHASE

2

## Surgical Solutions

- Voice Prosthesis Placement Surgical Kits
- Sterile Option Voice Prostheses
- Fistula Prostheses

PHASE

3

## Post-Operative Care

- Post-Operative Care Kits
- Laryngectomy Tubes
- Electrolarynx
- "Life After Laryngectomy" Brochure

PHASE

4

## Next Steps

- Next Steps Kit
- Voice Prostheses
- HMEs
- Adhesive Housings
- Educational Programs

PHASE

5


## And Beyond

- Hands-Free Speech Valve
- Skin Care

# ● PRE-OPERATIVE EDUCATION

## “What to Expect” Brochure

Available online and free to download at:  
<https://inhealth.com/educational-materials>

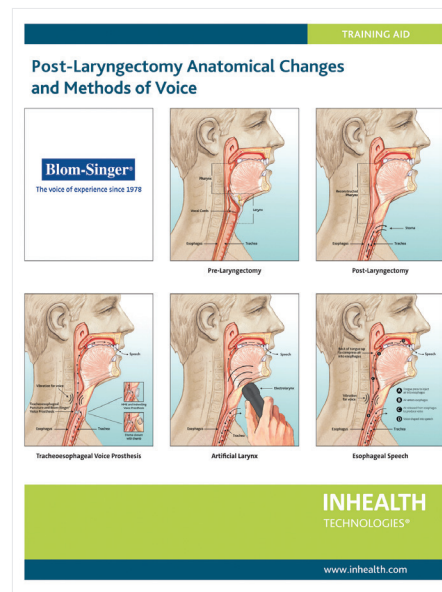

## Educational Videos

Available at InHealth Technologies®  
<https://inhealth.com/videos>



## ● SURGICAL SOLUTIONS

### The Blom-Singer® FirstFit™ Surgical Kit (Sterile)

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

Designed for primary or secondary tracheoesophageal puncture (TEP) with retrograde voice prosthesis insertion. This set includes: one (1) 13-gauge curved cored puncture needle; one (1) 36-inch guidewire with rounded tips; one (1) pharynx protector tool with handle and notched cylinder; one (1) voice prosthesis inserter body pre-loaded with a Blom-Singer Indwelling Voice Prosthesis of the diameter and length specified on the package label.

Ref	Description
TP 1608	FirstFit Surgical Kit + Classic 16 Fr 8 mm Voice Prosthesis
TP 1610	FirstFit Surgical Kit + Classic 16 Fr 10 mm Voice Prosthesis
TP 1612	FirstFit Surgical Kit + Classic 16 Fr 12 mm Voice Prosthesis
TP 1614	FirstFit Surgical Kit + Classic 16 Fr 14 mm Voice Prosthesis
TP 2008	FirstFit Surgical Kit + Classic 20 Fr 8 mm Voice Prosthesis
TP 2010	FirstFit Surgical Kit + Classic 20 Fr 10 mm Voice Prosthesis
TP 2012	FirstFit Surgical Kit + Classic 20 Fr 12 mm Voice Prosthesis
TP 2014	FirstFit Surgical Kit + Classic 20 Fr 14 mm Voice Prosthesis



**NEW**

### Blom-Singer® Tracheoesophageal Puncture Set (Sterile)

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

Designed for secondary placement of a tracheoesophageal prosthesis. Includes a disposable 16 gauge curved introducer needle with a 14 gauge Teflon® sheath, a 26" (66cm) silicone catheter with a tapered reinforced tip attached to a 24" (61cm) Teflon-coated stainless steel leader cable, and a fenestrated surgical drape with adhesive around the aperture.

Ref	Description
BE 6055	Tracheoesophageal Puncture Set (Sterile)



### Blom-Singer® Classic™ Indwelling Voice Prosthesis (Sterile)

Designed for primary TE puncture and placement procedures. Includes insertion and cleaning accessories. Clinician placed.

Length	16 Fr	20 Fr	Length	16 Fr	20 Fr
8mm	IN 1608-SO	IN 2008-SO	12mm	IN 1612-SO	IN 2012-SO
10mm	IN 1610-SO	IN 2010-SO	14mm	IN 1614-SO	IN 2014-SO




### Blom-Singer® Adjustable Bi-Flanged Fistula Prosthesis (Sterile)

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

Designed by Eric D. Blom, Ph.D to control leakage of saliva, food, and drink around a hypopharyngeal fistula. The prosthesis is adjustable in length and comes in two flange sizes, 25mm and 38mm, ensuring a proper fit.

Ref	Description
FM 1003	25mm (Small) Adjustable Bi-Flanged Fistula Prosthesis, 1 ea
FM 1004	38mm (Medium) Adjustable Bi-Flanged Fistula Prosthesis, 1 ea

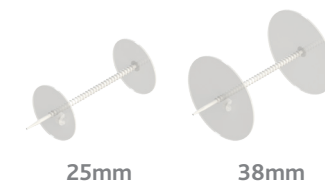

25mm

38mm

## ● POST-OPERATIVE CARE

### Blom-Singer® StomaSoft® Laryngectomy Tubes

Blom-Singer StomaSoft Laryngectomy Tubes are flexible, soft, and compatible with HME and other accessories from InHealth Technologies®. They feature reinforced strap holders and are available with or without fenestrations.

Size/Length	I.D./O.D.	Non-Fenestrated	Fenestrated	Rx
8/36	9.5mm/12.0mm	BE 6698	BE 6698F	■
8/55	9.5mm/12.0mm	BE 6699	BE 6699F	■
9/36	10.5mm/13.5mm	BE 6700	BE 6700F	■
9/55	10.5mm/13.5mm	BE 6701	BE 6701F	■
10/36	12.0mm/15.0mm	BE 6702	BE 6702F	■
10/55	12.0mm/15.0mm	BE 6703	BE 6703F	■
12/36	13.5mm/17.0mm	BE 6704	BE 6704F	■
12/55	13.5mm/17.0mm	BE 6705	BE 6705F	■



### Blom-Singer® Laryngectomy Tubes (Sterile)

Blom-Singer Laryngectomy Tubes (Sterile) are durable, lightweight silicone tubes that are used to keep the tracheostoma open after laryngectomy. The tubes do not hold an HME.

Size/Length	I.D./O.D.	Non-Fenestrated (Sterile)	Rx
8/36	9.5mm/12.0mm	BE 6398	■
8/55	9.5mm/12.0mm	BE 6399	■
9/36	10.5mm/13.5mm	BE 6400	■
9/55	10.5mm/13.5mm	BE 6401	■
10/36	12.0mm/15.0mm	BE 6402	■
10/55	12.0mm/15.0mm	BE 6403	■
12/36	13.5mm/17.0mm	BE 6404	■
12/55	13.5mm/17.0mm	BE 6405	■



### Marpac Tracheostomy Tube Holder/Collar™

Marpac Tracheostomy Tube Holder/Collar (for laryngectomy patients) is a single use, disposable Comfortable Collar™ used to secure the tube. Latex free. Adult—one size fits all.

Ref	Description
204	Tracheostomy Tube Holder/Collar



### StomaSoft Cleaning Swabs

The StomaSoft Cleaning Swabs are intended for cleaning tubes after removal from the patient. Available in packs of 10.

Ref	Description
TSWB10	Foam Swab, pack of 10



## Blom-Singer® Electrolarynx EL 1000

An electrolarynx is a speaking aid that produces sound vibrations that can be conducted through the throat and mouth to speak after a laryngectomy.

InHealth Technologies® and Eric D. Blom, Ph.D are pleased to offer a high quality, cost effective electrolarynx for laryngectomees.

Ref	Description
EL 1000	Electrolarynx Digital Speech Aid (One 9-volt batteries included)
EL 1002	Oral Adapter

Oral Adapter EL 1002



Electrolarynx EL 1000



Mark  
Year of TL: 2012  
Wearing AccuFit housing  
and Day&Night HME



## Blom-Singer® Post-Operative Care Kits (For Total Laryngectomy)

**CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN**

The Blom-Singer® Post-Operative Care Kits were designed with valued input from clinicians and patients. These kits contain supplies to assist patients in their transition from the hospital to home after a total laryngectomy. They are available in several configurations to meet select needs. All configurations include your choice of StomaSoft Laryngectomy Tube and HMEs as listed below.

The Post-Operative Care Kit is housed in a hard, clear storage case, which also includes information on how to obtain our Next Steps kit. The Next Steps kit includes additional sample product and educational items.

All kit configurations contain the following patient care supplies:

- 7 ea. 204 Marpac Tracheostomy Tube Holder/Collar™
- BE 6048 Blom-Singer® Shower Guard
- Drawstring Bag
- HydroFit Adhesive Housings, 6 pack
- Mirror with Stand
- Halcyon® LED Flashlight
- Next Steps Kit - Request Flyer

REF	Description
<b>PK 2100</b>	
BE 6701	Blom-Singer® StomaSoft® Laryngectomy Tube 9/55 Non-fenestrated, 2 ea
BE 1091	Day&Night HME® Cartridge ClassicFlow Heat and Moisture Exchange, 1 ea
<b>PK 2101</b>	
BE 6703	Blom-Singer® StomaSoft® Laryngectomy Tube 10/55 Non-fenestrated, 2 ea
BE 1091	Day&Night HME® Cartridge ClassicFlow Heat and Moisture Exchange, 1 ea
<b>PK 2102</b>	
BE 6701	StomaSoft Laryngectomy Tube 9/55 Non-fenestrated, 2 ea
BE 1092	Day&Night HME® Cartridge EasyFlow Heat and Moisture Exchange, 1 ea
<b>PK 2103</b>	
BE 6703	Blom-Singer® StomaSoft® Laryngectomy Tube 10/55 Non-fenestrated, 2 ea
BE 1092	Day&Night HME® Cartridge EasyFlow Heat and Moisture Exchange, 1 ea
<b>PK 2104</b>	
EL 1000	Blom-Singer® Electrolarynx, 1 ea
BE 6703	Blom-Singer® StomaSoft® Laryngectomy Tube 10/55 Non-fenestrated, 2 ea
BE 1091	Day&Night HME® Cartridge ClassicFlow Heat and Moisture Exchange, 1 ea
<b>PK 2106</b>	
BE 6705	Blom-Singer® StomaSoft® Laryngectomy Tube 12/55 Non-fenestrated, 2 ea
BE 1091	Day&Night HME® Cartridge ClassicFlow Heat and Moisture Exchange, 1 ea
<b>PK 2107</b>	
BE 6705	Blom-Singer® StomaSoft® Laryngectomy Tube 12/55 Non-fenestrated, 2 ea
BE 1092	Day&Night HME® Cartridge EasyFlow Heat and Moisture Exchange, 1 ea





Mike  
Year of TL: 2015  
Wearing StomaSoft  
laryngectomy tube and  
Day&Night HME



## ● NEXT STEPS

### VOICE PROSTHESES (CLINICIAN PLACED)

#### Blom-Singer® Classic™ Indwelling Voice Prosthesis (Sterile)

Includes sterile prostheses for primary TE puncture and placement procedures. Includes insertion and cleaning accessories. Clinician placed.

Length	16 Fr	20 Fr	Rx	Length	16 Fr	20 Fr	Rx
8mm	IN 1608-SO	IN 2008-SO	■	12mm	IN 1612-SO	IN 2012-SO	■
10mm	IN 1610-SO	IN 2010-SO	■	14mm	IN 1614-SO	IN 2014-SO	■



#### Blom-Singer® Classic™ Indwelling Voice Prosthesis

Includes nonsterile voice prosthesis for secondary or replacement TEP procedures, with insertion/cleaning accessories. Clinician placed.


Length	16 Fr	20 Fr	Rx	Length	16 Fr	20 Fr	Rx
4mm	IN 1604-NS	IN 2004-NS	■	14mm	IN 1614-NS	IN 2014-NS	■
6mm	IN 1606-NS	IN 2006-NS	■	16mm	IN 1616-NS	IN 2016-NS	■
8mm	IN 1608-NS	IN 2008-NS	■	18mm	IN 1618-NS	IN 2018-NS	■
10mm	IN 1610-NS	IN 2010-NS	■	20mm	IN 1620-NS	IN 2020-NS	■
12mm	IN 1612-NS	IN 2012-NS	■				



#### Blom-Singer® Dual Valve™ Indwelling Voice Prosthesis

The Blom-Singer Dual Valve Indwelling Voice Prosthesis has two valves embedded with silver oxide as a material preservative. When the primary valve fails, the secondary speech valve is intended to remain viable to help prevent leakage of fluids. The Blom-Singer® Dual Valve Voice Prosthesis is available with standard or large tracheal/esophageal flanges. Includes nonsterile voice prosthesis for secondary or replacement TEP procedures, with insertion/cleaning accessories. Clinician placed.

Length	20 Fr	20 Fr Large Flange	Rx	Length	20 Fr	20 Fr Large Flange	Rx
6mm	DV 2006-NS	DV 2006-LF	■	12mm	DV 2012-NS	DV 2012-LF	■
8mm	DV 2008-NS	DV 2008-LF	■	14mm	DV 2014-NS	DV 2014-LF	■
10mm	DV 2010-NS	DV 2010-LF	■				

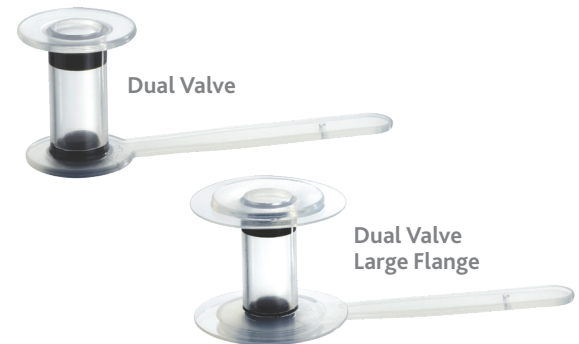

#### Blom-Singer® Advantage® Indwelling Voice Prosthesis


The Blom-Singer Advantage Indwelling Voice Prosthesis is embedded with silver oxide as a material preservative.\* Includes nonsterile voice prosthesis for secondary or replacement TEP procedures, with insertion/cleaning accessories.

##### HARD VALVE ASSEMBLY

Length	20 Fr	Rx
4mm	AD 2004-NS	■
6mm	AD 2006-NS	■
8mm	AD 2008-NS	■
10mm	AD 2010-NS	■
12mm	AD 2012-NS	■
14mm	AD 2014-NS	■

##### SOFT VALVE ASSEMBLY

Length	16 Fr	20 Fr	Rx
6mm	AS 1606-NS	AS 2006-NS	■
8mm	AS 1608-NS	AS 2008-NS	■
10mm	AS 1610-NS	AS 2010-NS	■
12mm	AS 1612-NS	AS 2012-NS	■
14mm	AS 1614-NS	AS 2014-NS	■

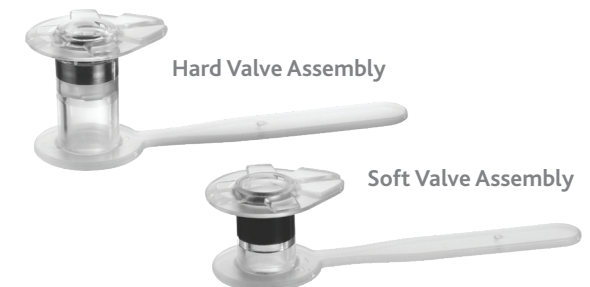

\*Silver oxide, a known anti-microbial substance, has been incorporated into the silicone flap valve of the Blom-Singer Advantage and Dual Valve Indwelling Voice Prostheses as a material preservative.



Left: Janice  
Year of TL: 2014  
Wearing AccuFit  
housing with  
SpeakFree HME  
Hands-Free Valve

Right: Mackayla  
Year of TL: 2021  
Wearing HydroFit  
housing and  
Day&Night HME

## SPECIAL ORDER VOICE PROSTHESES (CLINICIAN PLACED)


—DEVELOPED BY ERIC D. BLOM, PH.D

Unique problems call for unique solutions. InHealth Technologies® offers several voice prostheses to address those needs.

### Special Length

Available in one millimeter incremental lengths for a more precise TE puncture fit.


	5mm	7mm	9mm	Rx
16 Fr	IN 1605-SL	IN 1607-SL	IN 1609-SL	■
20 Fr	IN 2005-SL	IN 2007-SL	IN 2009-SL	■

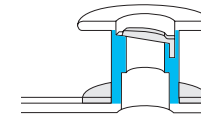

### Increased Resistance

Designed for patients who experience intermittent valve opening with fluctuations in esophageal pressures, resulting in air ingestion and/or leakage.


	6mm	8mm	10mm	12mm	Rx
16 Fr	IN 1606-IR	IN 1608-IR	IN 1610-IR	IN 1612-IR	■
20 Fr	IN 2006-IR	IN 2008-IR	IN 2010-IR	IN 2012-IR	■

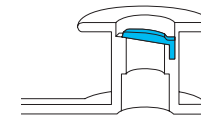

### Large Esophageal Flange


Provides a seal to help manage leakage around a voice prosthesis. Also, it enables downsizing from a wider diameter prosthesis to a smaller diameter. Effectiveness depends on an optimal prosthesis length that achieves a snug flange-to-flange seal against the wall surface(s).

	4mm	6mm	8mm	10mm	12mm	14mm	Rx
16 Fr	IN 1604-LEF	IN 1606-LEF	IN 1608-LEF	IN 1610-LEF	IN 1612-LEF	IN 1614-LEF	■
20 Fr	IN 2004-LEF	IN 2006-LEF	IN 2008-LEF	IN 2010-LEF	IN 2012-LEF	IN 2014-LEF	■

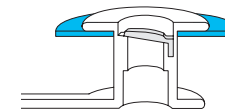

### LARGE ESOPHAGEAL FLANGE (SPECIAL LENGTH)

	5mm	7mm	Rx
16 Fr	IN 1605-LESL	IN 1607-LESL	■
20 Fr	IN 2005-LESL	IN 2007-LESL	■

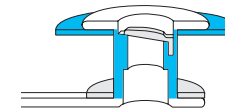

## Large Esophageal and Tracheal Flanges

Provides a double seal in patients with enlarged puncture tracts and tear-drop or irregularly shaped punctures. Effectiveness depends on an optimal prosthesis length fit that achieves a snug flange-to-flange seal against the wall surface(s).

	4mm	6mm	8mm	10mm	12mm	14mm	Rx
16 Fr	IN 1604-LF	IN 1606-LF	IN 1608-LF	IN 1610-LF	IN 1612-LF	IN 1614-LF	■
20 Fr	IN 2004-LF	IN 2006-LF	IN 2008-LF	IN 2010-LF	IN 2012-LF	IN 2014-LF	■




### LARGE ESOPHAGEAL FLANGE (INCREASED RESISTANCE)

	4mm	6mm	8mm	Rx
16 Fr	IN 1604-LEIR	IN 1606-LEIR	IN 1608-LEIR	■
20 Fr	IN 2004-LEIR	IN 2006-LEIR	IN 2008-LEIR	■

## TEP Occluders

The Blom-Singer® TEP (Tracheoesophageal Puncture) Occluder is a non-valved silicone device designed for placement in, and maintenance of, the TEP following total laryngectomy. Not a voice prosthesis.

	4mm	6mm	8mm	10mm	12mm	14mm	Rx
16 Fr	IN 1604-LETO	IN 1606-LETO	IN 1608-LETO	IN 1610-LETO	IN 1612-LETO	IN 1614-LETO	■
20 Fr	IN 2004-LETO	IN 2006-LETO	IN 2008-LETO	IN 2010-LETO	IN 2012-LETO	IN 2014-LETO	■

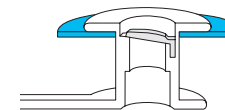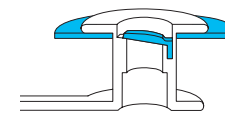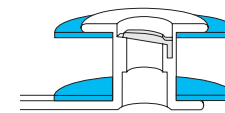


## VOICE PROSTHESES (PATIENT CHANGEABLE)

### Blom-Singer® Duckbill Voice Prostheses

The Blom-Singer® Duckbill is a patient changeable, slit valve prosthesis.


Length	16 Fr	Rx	Length	16 Fr	Rx
6mm	DB 16-006	■	12mm	DB 16-012	■
8mm	DB 16-008	■	14mm	BE 6001	■
10mm	BE 6000	■	18mm	BE 6002	■

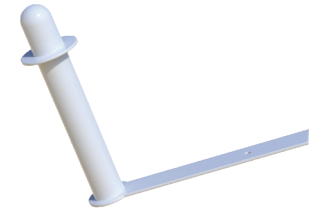

### Blom-Singer® Low Pressure Voice Prostheses

The Blom-Singer® Low Pressure is a patient changeable, one-way valve prosthesis.

Length	16 Fr	20 Fr	Rx	Length	16 Fr	20 Fr	Rx
6mm	BE 6009	BE 2014	■	18mm	BE 6012	BE 2026	■
8mm	LP 16-008	LP 20-008	■	22mm	BE 6013	BE 2030	■
10mm	BE 6010	BE 2018	■	25mm	BE 6014	BE 2033	■
12mm	LP 16-012	LP 20-012	■	28mm	BE 6015	—	■
14mm	BE 6011	BE 2022	■				

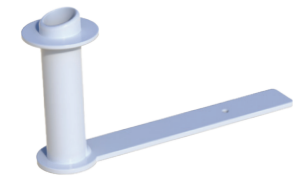

### Blom-Singer® Special Order Low Pressure Voice Prostheses

The Blom-Singer® Special Order Low Pressure is a patient changeable, one-way valve prosthesis.

Length	16 Fr Special Increased Resistance	16 Fr Special Length	20 Fr Special Length	20 Fr Special Increased Resistance	Rx
4mm	—	SPECIAL 11-902-06	SPECIAL 11-251-06	—	■
6mm	LP 16-006-IR	—	—	LP 20-006-IR	■
8mm	LP 16-008-IR	—	—	LP 20-008-IR	■
10mm	LP 16-010-IR	—	—	LP 20-010-IR	■
12mm	LP 16-012-IR	—	—	LP 20-012-IR	■
14mm	LP 16-014-IR	—	—	LP 20-014-IR	■
16mm	—	SPECIAL 11-902-04	SPECIAL 11-251-04	—	■
18mm	LP 16-018-IR	—	—	—	■
20mm	—	SPECIAL 11-902-05	SPECIAL 11-251-05	—	■
22mm	LP 16-022-IR	—	—	—	■
25mm	LP 16-025-IR	—	—	—	■
28mm	LP 16-028-IR	—	—	—	■

## Blom-Singer® Voice Prosthesis Sizers

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

The Blom-Singer Voice Prosthesis Sizer is a device used by clinicians to measure the length of the tracheoesophageal puncture to assist in selection of the appropriate voice prosthesis.

Ref	Description
BE 6036	Voice Prosthesis Sizer with 3 inserters (for 16 Fr voice prostheses)
BE 2036	Voice Prosthesis Sizer with 3 inserters (for 20 Fr voice prostheses)
BE 6096	Replacement Inserters (for use with both BE 6036 and BE 2036 Sizers), 10 ea



## Blom-Singer® Dilator-Sizers

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

The Blom-Singer Dilator-Sizer combines two tasks into one innovative product: sizing and dilation. In one step, the clinician can dilate and measure the tracheoesophageal puncture during voice prosthesis placement procedures.

Ref	Description
BE 2051	18 Fr Dilator-Sizer (for 16 Fr prostheses)
BE 2052	22 Fr Dilator-Sizer (for 20 Fr prostheses)



## VOICE PROSTHESIS ACCESSORIES

### Blom-Singer® Replacement Gel Caps

Use of the Blom-Singer Replacement Gel Caps provides a smooth, rounded shape to the esophageal flange end of the voice prosthesis, enabling easy insertion into the tracheoesophageal puncture.

Ref	Description	Rx
BE 3190	Replacement Gel Caps, 16 Fr, 90 ea	■
BE 3195	Replacement Gel Caps, 18 Fr, 90 ea	■
BE 3290	Replacement Gel Caps, 20 Fr, 90 ea	■
BE 3295	Replacement Gel Caps, 22 Fr, 60 ea	■



### Blom-Singer® Flushing Device and Cleaning Brushes

For routine cleaning of a voice prosthesis as recommended by your clinician.

Ref	Description
IN 4050	Flushing Device, pack of 3
IN 4064*	Cleaning Brush, 16/20 Fr, Small, 4mm–8mm, 3 ea
IN 4065*	Cleaning Brush, 16/20 Fr, Medium, 9mm–16mm, 3 ea
IN 4066*	Cleaning Brush, 16/20 Fr, Large, 18mm–22mm, 3 ea

\*Fits indwelling style voice prostheses only.



### Blom-Singer® Plug Inserts

The Blom-Singer Plug Insert is designed for temporary management of leakage through a voice prosthesis. Not compatible with patient changeable voice prostheses.

Ref	Description	Rx
IN 4069	16 Fr Plug Insert (Advantage and Classic, 6mm to 20mm)	■
IN 4070	20 Fr Plug Insert (Advantage and Classic, 6mm to 20mm)	■



# Blom-Singer® HMEs

## Simple choices to fit your lifestyle

Left: Joanne  
Year of TL: 2015  
Wearing AccuFit  
housing with  
SpeakFree HME  
Hands-Free

Right: Meagon  
Year of TL: 2016  
Wearing HydroFit  
housing with  
Day&Night HME

### Day&Night™ HME

- Wear all day and all night, including bedtime
- Attaches to the entire family of InHealth baseplates and StomaSoft tubes.



Day&Night®  
EasyFlow® HME and  
ClassicFlow® HME

### SpeakFree™ HME

- Your HME and hands free speaking valve, all in one.
- Attaches to the entire family of InHealth baseplates and StomaSoft tubes.



SpeakFree®  
EasyFlow® HME and  
ClassicFlow® HME





## HEAT & MOISTURE EXCHANGE (HME)

Research shows that long-term HME use results in less mucus<sup>1</sup>, less coughing<sup>2</sup>, and decreased shortness of breath, fatigue, sleep problems, anxiety, and depression<sup>1,2</sup>.

### Blom-Singer® SpeakFree® Hands-Free Valve with ClassicFlow® HME

Blom-Singer SpeakFree HME Hands-Free Valve is the first heat and moisture exchange (HME) cartridge with a single-use fully integrated adjustable hands-free valve. The SpeakFree HME is translucent and non-reflective for a natural look. The SpeakFree ClassicFlow HME is designed for increased moisture retention and resistance.

Ref	Description
BE 1090EZ	SpeakFree Hands-Free Valve with ClassicFlow HME, 30 per pack



SpeakFree  
ClassicFlow

### Blom-Singer® SpeakFree® Hands-Free Valve with EasyFlow® HME

Blom-Singer SpeakFree HME Hands-Free Valve is the first heat and moisture exchange (HME) cartridge with a single-use fully integrated adjustable hands-free valve. The SpeakFree HME is translucent and non-reflective for a natural look. With SpeakFree EasyFlow, the foam filter is more porous and allows you to breathe more freely to match your activity level and pulmonary needs.

Ref	Description
BE 1090EF	SpeakFree Hands-Free Valve with EasyFlow HME, 30 per pack



SpeakFree  
EasyFlow

### Blom-Singer® Day&Night® ClassicFlow® HME

The Blom-Singer Day&Night ClassicFlow can be worn all day and night including bedtime. The Day&Night HME is translucent and non-reflective for a natural look.

Ref	Description
BE 1091	Day&Night ClassicFlow HME Cartridge, pack of 30 (24-hour use)



Day&Night  
ClassicFlow

### Blom-Singer® Day&Night® EasyFlow® HME

The Blom-Singer Day&Night EasyFlow allows you to breathe more freely to match your activity level and pulmonary needs. The Day&Night HME is translucent and non-reflective for a natural look.

Ref	Description
BE 1092	Day&Night EasyFlow HME Cartridge, pack of 30 (24-hour use)



Day&Night  
EasyFlow

1. Ackerstaff AH, Hilgers FJ, et al. Improvements in respiratory and psychosocial functioning following total laryngectomy by the use of a heat and moisture exchanger. *Annals of Otolaryngology, Rhinology & Laryngology*. 1993;102(11):878-883.

2. Ackerstaff AH, Hilgers FJ, et al. Heat and moisture exchangers as a treatment option in the post-operative rehabilitation of laryngectomized patients. *Clinical Otolaryngology*. 1995;20:504-509.

## ADHESIVE HOUSINGS

InHealth Technologies® presents the FitSeries—adhesive housings designed to meet your skin and stoma needs after a total laryngectomy. Available in two styles, three shapes each. Compatible with Day&Night and SpeakFree Hands Free HME cartridges.

### Blom-Singer® AccuFit® Adhesive Housings

The Blom-Singer AccuFit Adhesive Housing is designed to be our strongest adhesive option while still maintaining flexibility. It has increased elasticity and is less prone to tearing. Available in three shapes: round, oval, and oval extra.


Ref	Description
BE 6082	AccuFit Adhesive Housing, Round, pack of 30
BE 6083	AccuFit Adhesive Housing, Oval, pack of 30
BE 6084	AccuFit Adhesive Housing, Oval Extra, pack of 30



### Blom-Singer® HydroFit® Adhesive Housings

The Blom-Singer HydroFit Adhesive Housing is a hydrocolloid adhesive baseplate made with a protective gel-like layer and a moisture-absorbent adhesive. The HydroFit is gentle on the skin and flexible, but maintains strong adhesion. Available in three shapes: round, oval, and oval extra.

Ref	Description
BE 6085	Hydrofit Adhesive Housing, Round, pack of 30
BE 6086	Hydrofit Adhesive Housing, Oval, pack of 30
BE 6087	Hydrofit Adhesive Housing, Oval Extra, pack of 30

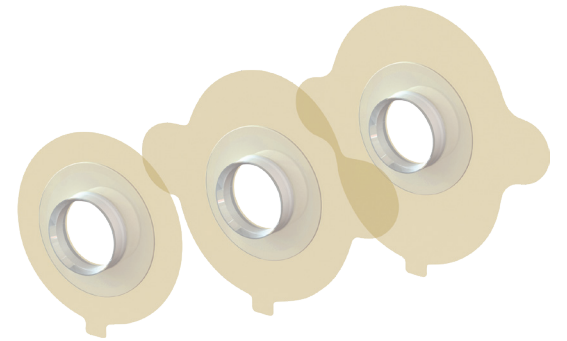

## TRACHEOSTOMA ATTACHMENTS

### Blom-Singer® Tracheostoma Valve Housings

The Tracheostoma Valve Housing is a reusable, flexible housing that attaches to the Day&Night and SpeakFree Hands-Free HME cartridge or Shower Guard over the stoma. It is applied with double sided tape discs. It is available in two sizes, standard and large.

Ref	Description
BE 6038	Tracheostoma Valve Housing, standard, 1 ea
BE 6039	Tracheostoma Valve Housing, large, 1 ea



### Blom-Singer® Adhesive Tape Discs

A variety of double-sided tape discs are available to provide attachment of the tracheostoma valve housing.

Ref	Description	Ref	Description
BE 6041	Regular Tape Discs, standard, pack of 30	BE 6044	Thin Tape Discs, large, pack of 30
BE 6042	Regular Tape Discs, large, pack of 30	BE 6034	Heavy Duty Tape Discs, standard, pack of 30
BE 6043	Thin Tape Discs, standard, pack of 30	BE 6035	Heavy Duty Tape Discs, large, pack of 30



Ksenia  
Year of TL: 2019  
Wearing Accufit  
housing and  
Day&Night HME



## ● AND BEYOND

### HANDS-FREE SPEECH VALVE

#### Blom-Singer® SpeakFree® Hands-Free Valve with ClassicFlow® HME

Blom-Singer SpeakFree HME Hands-Free Valve is the first heat and moisture exchange (HME) cartridge with a single-use fully integrated adjustable hands-free valve. The SpeakFree HME is translucent and non-reflective for a natural look. The SpeakFree ClassicFlow HME is designed for increased moisture retention and resistance.

Ref	Description
BE 1090EZ	SpeakFree Hands-Free Valve with ClassicFlow HME, 30 per pack

#### Blom-Singer® SpeakFree® Hands-Free Valve with EasyFlow® HME

Blom-Singer SpeakFree HME Hands-Free Valve is the first heat and moisture exchange (HME) cartridge with a single-use fully integrated adjustable hands-free valve. The SpeakFree HME is translucent and non-reflective for a natural look. With SpeakFree EasyFlow, the foam filter is more porous and allows you to breathe more freely to match your activity level and pulmonary needs.

Ref	Description
BE 1090EF	SpeakFree Hands-Free Valve with EasyFlow HME, 30 per pack



SpeakFree  
ClassicFlow



SpeakFree  
EasyFlow

### SKIN CARE

#### Skin-Prep™ Liquid Film-Forming Protective Barrier Wipes

Skin-Prep is a liquid film-forming dressing that, upon application to intact skin, forms a protective film to help reduce friction during removal of tapes and films.

Ref	Description
59420425	Skin-Prep Liquid Film-Forming Protective Barrier Wipes, box of 50

#### Skin Tac™ Adhesive Barrier Wipes

Skin Tac liquid adhesive barrier wipes are latex-free, hypo-allergenic, clear, non-rubber adhesive barrier wipe.

Ref	Description
MS407W*	Adhesive Barrier Wipes, box of 50

\*May not be available in all international markets.

#### NFS Adhesive (Silicone Adhesive)

PSA Pressure Sensitive Adhesive, silicone based adhesive for use with all silicone prosthetic devices.

Ref	Description
B-530	NFS Adhesive, 1 oz

#### Remove™ Adhesive Remover

Gently cleans hydrocolloid, acrylic and rubber-based adhesive residues from the skin without irritation, discomfort or trauma.

Ref	Description
59403125	Adhesive Remover, box of 50



## STOMA PROTECTION AND ACCESSORIES

### ADDvox® Stoma Cover Mock Turtlenecks (Only available in the US)

The ADDvox Stoma Cover Mock Turtleneck has inner layers of material which helps keep debris from entering the stoma. They slip under a shirt like a fashion dickey and have an outer layer of soft 100% cotton knit. Each cover is 9-1/2" long x 8-3/4" wide. Velcro fastener adjusts to fit necks of 12-1/2" to 22". Each color sold separately.

Ref	Description	Color	Quantity
BE 6227-R2	Stoma Cover Mock Turtleneck	White	1 ea
BE 6228-R2	Stoma Cover Mock Turtleneck	Black	1 ea
BE 6229-R2	Stoma Cover Mock Turtleneck	Navy	1 ea



### ADDvox® Stoma Scarves (Only available in the US)

The ADDvox Stoma Cover Scarf has inside layers to help filter out dust, dirt and pollen, is completely washable, and needs little care. This attractive neckline scarf looks just like a regular scarf or ascot and is available in a wide range of colors for men or women. All scarves have Velcro closures and adjust to fit necks from 17-3/4" to 22" (can be made smaller with tucks or stitching). Each color sold separately.

Ref	Description	Color	Quantity
BE 6232-R2	Stoma Cover Scarf	Black	1 ea
BE 6233-R2	Stoma Cover Scarf	White	1 ea
BE 6234-R2	Stoma Cover Scarf	Navy	1 ea



### ADDvox® Cloth Bib-Style Stoma Covers (Only available in the US)

The ADDvox soft, long-lasting bib protects the tracheostoma and lets you breathe more easily. This "bib-style" Stoma Cover will look crisp and clean after repeated laundering. Five soft layers keep dust and dirt from reaching the trachea. Fits necks up to 22". Velcro closure. Each color sold separately.

Ref	Description	Color	Quantity
BE 6224-R2	Cloth Bib-Style Stoma Cover	White	1 ea
BE 6225-R2	Cloth Bib-Style Stoma Cover	Beige	1 ea
BE 6226-R2	Cloth Bib-Style Stoma Cover	Navy	1 ea



## ADDvox® Silicone Shower Shield (Only available in the US)

The one-piece ADDvox Silicone Shower Shield is worn during showering to deflect the water away from your stoma. Extremely soft and pliable, it creates a soft, comfortable and conforming seal along the irregular contours of your neck, and vents on the underside allowing for airflow. Made of medical grade silicone, this latex-free Silicone Shower Shield is comfortable, durable, and not affected by temperature, ozone, sunlight, soap or detergents.

Ref	Description
BE 6079-R2	Silicone Shower Shield, 1 ea



## Blom-Singer® Shower Guard

The Blom-Singer Shower Guard is a one-piece, right-angle plastic shower guard that comes with a Tracheostoma Standard Adhesive Housing and 30 double-sided Adhesive Tape Discs. Also compatible with the Blom-Singer StomaSoft® Laryngectomy Tube.

Ref	Description
BE 6048	Shower Guard—Tracheostoma Valve Housing (BE 6038) and Regular Tape Discs, pk 30 (BE 6041) included



## Blom-Singer® Foam Stoma Protector

A simple and temporary (disposable) solution to help prevent debris from entering the stoma. Adhesive strip on one edge allows for easy attachment, access to the stoma to occlude for speech with a voice prosthesis, or for coughing. Measurement listed is approximate.

Ref	Description
BE 6220	Foam Stoma Protector, Tan, Large 2.7" L x 2.7" W (7cm x 7cm), pack of 30



## ADDvox® Stoma Filters with Microporous Adhesive

ADDvox Stoma Filters offer protection against dust, dirt and pollen. Designed for single use with a hypoallergenic, non-irritating, self-adhesive microporous strip that is comfortable and easy to apply. Not recommended for cannula users. 2-1/4" wide x 2-1/2" high x 1/8" thick. Available in off-white and skin tone.

Ref	Description
BE 6222-R2	Stoma Filters, Skin tone, pack of 60
BE 6223-R2	Stoma Filters, Off-white, pack of 60





**Professional Products**

## VOICE PROSTHESIS MANAGEMENT

### Blom-Singer® Voice Prosthesis Sizers

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

The Blom-Singer Voice Prosthesis Sizer is a device used by clinicians to measure the length of the tracheoesophageal puncture to assist in selection of the appropriate voice prosthesis.

Ref	Description
BE 6036	Voice Prosthesis Sizer with 3 inserters (for 16 Fr voice prostheses)
BE 2036	Voice Prosthesis Sizer with 3 inserters (for 20 Fr voice prostheses)
BE 6096	Replacement Inserters (for use with both BE 6036 and BE 2036 Sizers), 10 ea



### Blom-Singer® Dilator-Sizers

CAUTION: FEDERAL (USA) LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN

The Blom-Singer Dilator-Sizer combines two tasks into one innovative product: sizing and dilation. In one step, the clinician can dilate and measure the tracheoesophageal puncture during voice prosthesis placement procedures.

Ref	Description
BE 2051	18 Fr Dilator-Sizer (for 16 Fr prostheses)
BE 2052	22 Fr Dilator-Sizer (for 20 Fr prostheses)



### Blom-Singer® Replacement Gel Caps

Ref	Description	Rx
BE 3190	Replacement Gel Caps, 16 Fr, 90 ea	■
BE 3195	Replacement Gel Caps, 18 Fr, 90 ea	■
BE 3290	Replacement Gel Caps, 20 Fr, 90 ea	■
BE 3295	Replacement Gel Caps, 22 Fr, 60 ea	■



### Manometer with Manometer Adaptor

A tool to teach patients how to self-regulate intratracheal air pressure. Proper airflow control may allow for better voice quality, less fatigue, and longer duration of the tracheostoma adhesive housing seal.

Ref	Description
BE 8000	Manometer Kit includes 1 manometer, 1 manometer adaptor and 1 adhesive housing
BE 8020	Replacement Manometer Adaptor includes 5 manometer adaptors and 5 adhesive housings



### Blom-Singer® Insufflation Test Set

The Blom-Singer Insufflation Test Set is designed to aid the clinician in subjectively predicting the ability of the laryngectomized patient to produce tracheoesophageal voice.

Ref	Description
BE 6040	Insufflation Test Set, box of 6







# The Choice for Long- or Short-Term Augmentation

## RENÚ® VOICE AND RENÚ® GEL

RENÚ Voice is for long-term augmentation:

- A buffered hydrogel of synthetic, proven-biocompatible calcium hydroxylapatite (CaHA) particles 25 to 45 microns in size for consistent and easy injection.<sup>1</sup>

RENÚ Gel is for short-term augmentation:

- A buffered hydrogel of sodium carboxymethylcellulose<sup>2</sup> (NaCMC) and glycerin.
- Resorbs slowly in vivo, providing a scaffold for local tissue infiltration.<sup>3</sup>

### RENÚ® Voice: Calcium Hydroxylapatite Implant

RENÚ Voice is an injectable vocal fold implant intended for long-term augmentation.

Ref	Description
08-015-04-V04	RENÚ Voice Calcium Hydroxylapatite Implant, 1.5cc
08-015-04-V01	RENÚ Voice Calcium Hydroxylapatite Implant, 1.5cc, CE

### RENÚ® Gel: Hydrogel Implant

RENÚ Gel is an injectable vocal fold implant intended for short-term augmentation.

Ref	Description
08-015-00-V04	RENÚ Gel Hydrogel Implant, 1.5cc

1. Data on file, InHealth Technologies®.

2. Mallur PS, et al. Safety and efficacy of carboxymethylcellulose in the treatment of glottic insufficiency. The Laryngoscope. 2012;122:322–326. <https://doi.org/10.1002/lary.21930>.

3. Hughes RGM, Morrison M. Vocal cord medialization by transcutaneous injection of calcium hydroxylapatite. Journal of Voice. 2005;19(4):674–678.

## RENÚ® VOICE



## RENÚ® GEL






**TO LEARN MORE**  
Scan QR code  
to visit the  
RENU website

# Katz Extractor<sup>®</sup>

OTO-RHINO FOREIGN BODY REMOVER

The easy way out

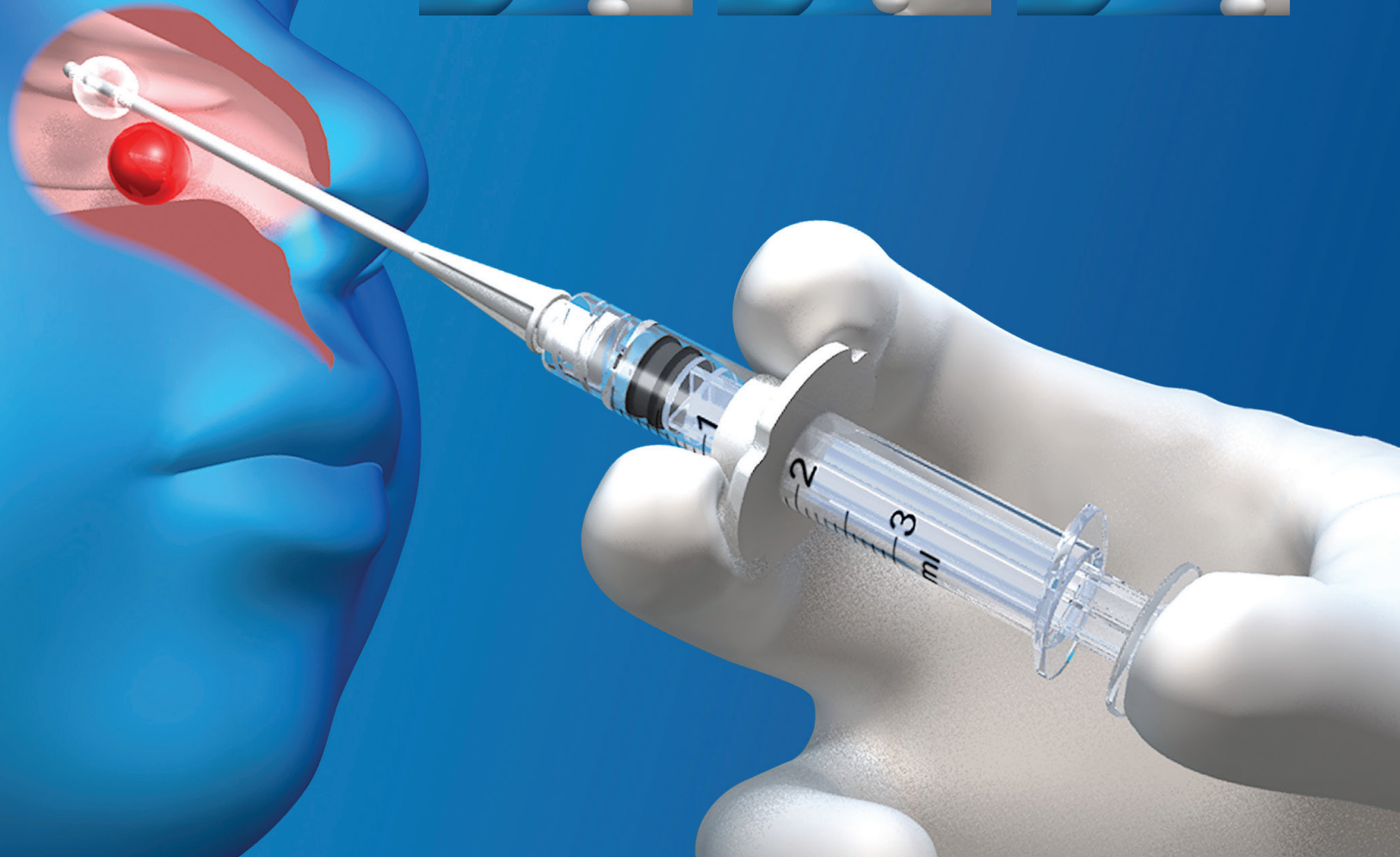


## Katz Extractor® Oto-Rhino Foreign Body Remover

Designed in conjunction with David L. Katz, MD, MPH, the Katz Extractor Oto-Rhino Foreign Body Remover is intended for extraction of impacted foreign bodies in the nasal passage. Single use only. This device may also be used in the external ear canal, based upon clinical judgment. Covered by one or more of the following patents and other patents pending: US#5,454,817. Katz Extractor® is a registered trademark of Freudenberg Medical, LLC.

Ref	Description
FB 5000	Oto-Rhino Foreign Body Remover, box of 5

**TO LEARN MORE**  
Scan QR code  
to visit the  
Katz Extractor website



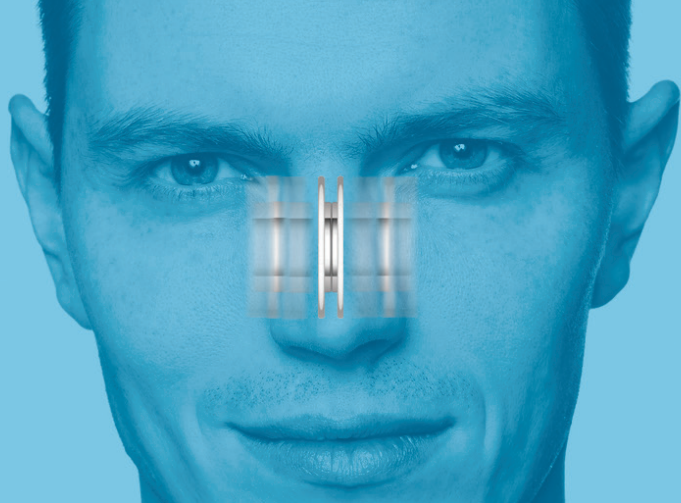

# Nasal Septal Perforation Prosthesis


A two-piece, magnetically coupled solution for non-surgical closure of nasal septal perforations.

## Nasal Septal Perforation Prostheses

The Blom-Singer® Nasal Septal Perforation Prosthesis is a two-piece, magnetically coupled solution for non-surgical closure of nasal septal perforations.








### OVAL

Ref	Description
SP 2112	Nasal Septal Perforation Prosthesis Oval 11 x 13 x 2mm
SP 2113	Nasal Septal Perforation Prosthesis Oval 11 x 13 x 3mm
SP 2132	Nasal Septal Perforation Prosthesis Oval 13 x 17 x 2mm
SP 2133	Nasal Septal Perforation Prosthesis Oval 13 x 17 x 3mm
SP 2152	Nasal Septal Perforation Prosthesis Oval 15 x 21 x 2mm
SP 2153	Nasal Septal Perforation Prosthesis Oval 15 x 21 x 3mm

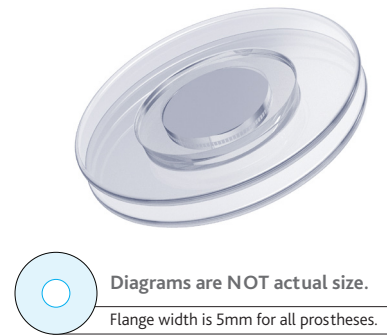

### ROUND

Ref	Description
SP 1052	Nasal Septal Perforation Prosthesis Round 5 x 2mm
SP 1053	Nasal Septal Perforation Prosthesis Round 5 x 3mm
SP 1072	Nasal Septal Perforation Prosthesis Round 7 x 2mm
SP 1073	Nasal Septal Perforation Prosthesis Round 7 x 3mm
SP 1102	Nasal Septal Perforation Prosthesis Round 10 x 2mm
SP 1103	Nasal Septal Perforation Prosthesis Round 10 x 3mm

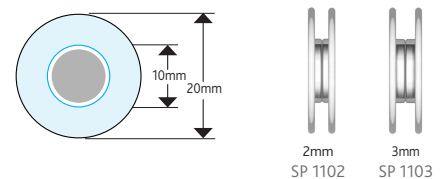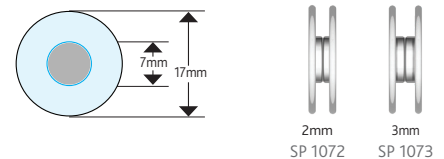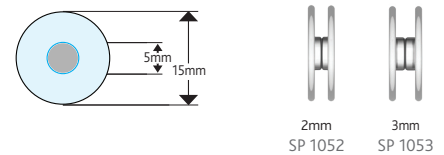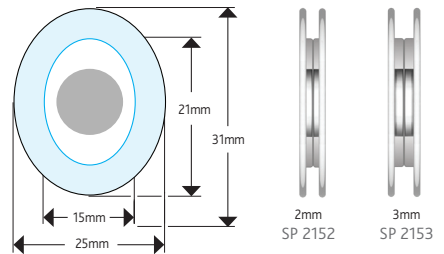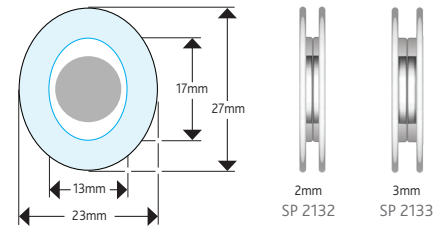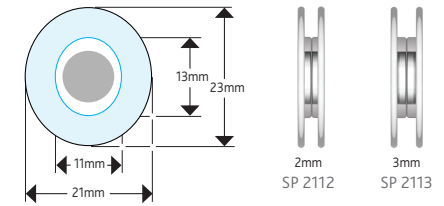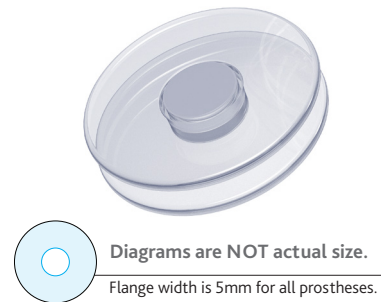

**TO LEARN MORE**  
Scan QR code to visit the  
Nasal Septal Perforation Prostheses website

## ORDERING INFORMATION

### UNITED STATES

**Phone:** Toll-Free: 800.477.5969

**Fax:** Toll-Free: 888.371.1530

**E-mail:** [order@inhealth.com](mailto:order@inhealth.com)

**Prescription E-mail:** [ihtrx@inhealth.com](mailto:ihtrx@inhealth.com)

**Order Online:** [www.inhealth.com](http://www.inhealth.com)

**Mail:**

InHealth Technologies®  
1110 Mark Ave.  
Carpinteria, CA 93013 • USA

### CANADA

**Phone:** Toll-Free: 1.888.437.4123

**Fax:** Toll-Free: 1.226.641.5285

**E-mail:** [canadaorder@inhealth.com](mailto:canadaorder@inhealth.com)

**Website:** [www.inhealth.com/canada](http://www.inhealth.com/canada)

**Mail:**

InHealth Canada  
65 Spruce Street  
Tillsonburg ON N4G 5C4

**Rx: Products  
with the Rx  
notation require  
a prescription  
from a physician  
(US patients).**

## NEW PATIENT ENROLLMENT INSTRUCTIONS

### Enrolling as a new patient? It's easy! Here is what to do:

1. Obtain a **prescription** from your provider.
2. Fill out the **Patient Enrollment Form**.
3. Send back to InHealth Technologies® with a copy of your insurance card (front and back). You can send via email: [medicare@inhealth.com](mailto:medicare@inhealth.com) or fax: **888.371.1530**
4. Place your order! Please remember that a prescription alone will not automatically generate an order.

### Need insurance or Medicare help? Contact us!

Phone: **800.477.5969 option 3**

Email: [medicare@inhealth.com](mailto:medicare@inhealth.com)

While you don't have to fill anything out on the **Patient Bill of Rights** form, please make sure you review it. We want you to know your rights.

All forms are available under the **Resources** tab on the InHealth Technologies® website: [www.inhealth.com](http://www.inhealth.com)

## FREE EDUCATIONAL RESOURCES

- Available on [www.inhealth.com](http://www.inhealth.com) under the **Resources** tab.
- Visit InHealth Technologies® video channel at <https://inhealth.com/videos> to view our videos

In 1978, Drs. Singer and Blom pioneered the methods of tracheoesophageal puncture and valved silicone voice prostheses that have been the international standard for voice restoration for over 45 years, allowing thousands of patients worldwide to regain their ability to speak. As a leading supplier of advanced voice restoration systems, our goal at InHealth Technologies® is to help laryngectomees regain and preserve a vital part of their lives—their voice. We are proud to continue the Blom-Singer® legacy today with clinical support, input, and leadership.


**Blom-Singer®**



Janice  
Year of TL: 2014  
Wearing AccuFit  
housing with  
SpeakFree HME  
Hands-Free Valve



## STAY CONNECTED...

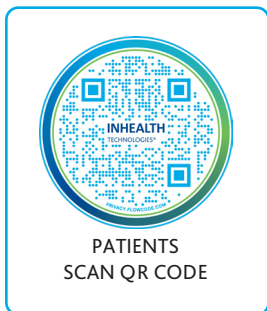

### Speakers Club™

**Patients: Scan the QR code to sign up today!**

We created this digital E-News for you! **Speakers Club E-News** will provide educational laryngectomy and voice prosthesis information.



### ProNET™

**Professionals: Scan the QR code to sign up today!**

We created this digital E-News specifically for professionals. **ProNET E-News** will deliver useful laryngectomy and voice prosthesis information for you and your practice, including material you may want to share with your patients.

**INHEALTH**  
TECHNOLOGIES®

**United States**  
1110 Mark Ave., Carpinteria, CA 93013  
Toll-Free: 800.477.5969  
Toll-Free Fax: 888.371.1530

**International**  
Tel: +1.805.684.9337  
Fax: +1.805.684.8594  
E-mail: [international@inhealth.com](mailto:international@inhealth.com)

Patents: <https://inhealth.com/media/documents/pdf/Patents.pdf>  
Terms & Conditions: [https://inhealth.com/media/documents/pdf/Terms\\_and\\_Conditions.pdf](https://inhealth.com/media/documents/pdf/Terms_and_Conditions.pdf)  
©2024 InHealth Technologies®—Manufactured by Freudenberg Medical, LLC (222002.02)

**Blom-Singer®**  
The voice of experience since 1978



Order online at [www.inhealth.com](http://www.inhealth.com)